

Archived Document

This archived document is no longer being reviewed through the CLSI Consensus Document Development Process. However, this document is technically valid as of January 2017. Because of its value to the laboratory community, it is being retained in CLSI's library.



October 2004

LIS02-A2

Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems; Approved Standard— Second Edition



This document covers the two-way digital transmission of remote requests and results between clinical laboratory instruments and information systems.

A standard for global application developed through the Clinical and Laboratory Standards Institute consensus process.

Clinical and Laboratory Standards Institute

Setting the standard for quality in medical laboratory testing around the world.

The Clinical and Laboratory Standards Institute (CLSI) is a not-for-profit membership organization that brings together the varied perspectives and expertise of the worldwide laboratory community for the advancement of a common cause: to foster excellence in laboratory medicine by developing and implementing medical laboratory standards and guidelines that help laboratories fulfill their responsibilities with efficiency, effectiveness, and global applicability.

Consensus Process

Consensus—the substantial agreement by materially affected, competent, and interested parties—is core to the development of all CLSI documents. It does not always connote unanimous agreement, but does mean that the participants in the development of a consensus document have considered and resolved all relevant objections and accept the resulting agreement.

Commenting on Documents

CLSI documents undergo periodic evaluation and modification to keep pace with advancements in technologies, procedures, methods, and protocols affecting the laboratory or health care.

CLSI's consensus process depends on experts who volunteer to serve as contributing authors and/or as participants in the reviewing and commenting process. At the end of each comment period, the committee that developed the document is obligated to review all comments, respond in writing to all substantive comments, and revise the draft document as appropriate.

Comments on published CLSI documents are equally essential, and may be submitted by anyone, at any time, on any document. All comments are managed according to the consensus process by a committee of experts.

Appeals Process

When it is believed that an objection has not been adequately considered and responded to, the process for appeals, documented in the CLSI Standards Development Policies and Processes, is followed.

All comments and responses submitted on draft and published documents are retained on file at CLSI and are available upon request.

Get Involved—Volunteer!

Do you use CLSI documents in your workplace? Do you see room for improvement? Would you like to get involved in the revision process? Or maybe you see a need to develop a new document for an emerging technology? CLSI wants to hear from you. We are always looking for volunteers. By donating your time and talents to improve the standards that affect your own work, you will play an active role in improving public health across the globe.

For additional information on committee participation or to submit comments, contact CLSI.

Clinical and Laboratory Standards Institute
950 West Valley Road, Suite 2500
Wayne, PA 19087 USA
P: +1.610.688.0100
F: +1.610.688.0700
www.clsi.org
standard@clsi.org

ISBN 1-56238-550-X
ISSN 0273-3099

LIS02-A2
Vol. 24 No. 33
Replaces LIS2-A
Vol. 23 No. 8

Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems; Approved Standard—Second Edition

Volume 24 Number 33

Paul J. Mountain, M.Sc., M.T.(ASCP)
David Chou, M.D.
James V. Callaghan, M.T.(ASCP)
Randall R. Davis
Charles D. Hawker, Ph.D., MBA, FACB
David A. Herold, M.D., Ph.D.
Andrzej J. Knafel, Ph.D.
Gary W. Kramer, Ph.D.
Rodney S. Markin, M.D., Ph.D.

Abstract

Clinical and Laboratory Standards Institute document LIS02-A2—*Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems; Approved Standard—Second Edition* addresses the two-way digital transmission of remote requests and results between clinical laboratory instruments and information systems. It enables any two such systems to establish a logical link for communicating text to send result, request, or demographic information in a standardized and interpretable form.

Clinical and Laboratory Standards Institute (CLSI). *Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems; Approved Standard—Second Edition*. CLSI document LIS02-A2 (ISBN 1-56238-550-X). Clinical and Laboratory Standards Institute, 950 West Valley Road, Suite 2500, Wayne, Pennsylvania 19087 USA, 2004.

The Clinical and Laboratory Standards Institute consensus process, which is the mechanism for moving a document through two or more levels of review by the health care community, is an ongoing process. Users should expect revised editions of any given document. Because rapid changes in technology may affect the procedures, methods, and protocols in a standard or guideline, users should replace outdated editions with the current editions of CLSI documents. Current editions are listed in the CLSI catalog and posted on our website at www.clsi.org. If your organization is not a member and would like to become one, and to request a copy of the catalog, contact us at: Telephone: 610.688.0100; Fax: 610.688.0700; E-Mail: customerservice@cls.org; Website: www.clsi.org.



Copyright ©2004 Clinical and Laboratory Standards Institute. Except as stated below, any reproduction of content from a CLSI copyrighted standard, guideline, companion product, or other material requires express written consent from CLSI. All rights reserved. Interested parties may send permission requests to permissions@clsi.org.

CLSI hereby grants permission to each individual member or purchaser to make a single reproduction of this publication for use in its laboratory procedure manual at a single site. To request permission to use this publication in any other manner, e-mail permissions@clsi.org.

Suggested Citation

CLSI. *Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems; Approved Standard—Second Edition*. CLSI document LIS02-A2. Wayne, PA: Clinical and Laboratory Standards Institute; 2004.

Previous Edition:

April 2003

Reaffirmed:

September 2011

Archived:

January 2017

ISBN 1-56238-550-X

ISSN 0273-3099

Committee Membership

Area Committee on Automation and Informatics

**Paul J. Mountain, M.Sc.,
M.T.(ASCP)
Chairholder
MDS Laboratories
Toronto, Ontario, Canada**

**David Chou, M.D.
Vice-Chairholder
University of Washington Medical
Center
Seattle, Washington**

James V. Callaghan, M.T.(ASCP)
FDA Center for Devices and
Radiological Health
Rockville, Maryland

Randall R. Davis
Dade Behring Inc.
Newark, Delaware

Charles D. Hawker, Ph.D., MBA,
FACB
ARUP Laboratories, Inc.
Salt Lake City, Utah

David A. Herold, M.D., Ph.D.
VA (San Diego) Medical Center
San Diego, California

Andrzej J. Knafel, Ph.D.
Roche Instrument Center AG
Rotkreuz, Switzerland

Gary W. Kramer, Ph.D.
National Institute of Standards &
Technology
Gaithersburg, Maryland

Rodney S. Markin, M.D., Ph.D.
Univ. of Nebraska Medical Center
Omaha, Nebraska

Advisors

Michael G. Bissell, M.D., Ph.D.,
M.P.H.
Ohio State University
Columbus, Ohio

Mary F. Burritt, Ph.D.
Mayo Clinic
Rochester, Minnesota

Suzanne H. Butch, M.A.,
M.T.(ASCP), SBB
The University of Michigan
Ann Arbor, Michigan

Al DeStefano
Sysmex Corporation
Tucson, Arizona

Robert J. Dominici
Cholestech Corporation
Alamo, California

Jeffrey A. DuBois, Ph.D.
Nova Biomedical Corporation
Waltham, Massachusetts

Louis J. Dunka, Jr., Ph.D.
LifeScan, Inc.
Milpitas, California

Robert H. Engel, Ph.D.
Engel Associates
Duxbury, Massachusetts

Arden W. Forrey, Jr., Ph.D., FACB
University of Washington
Seattle, Washington

Masayoshi Hayashi
Sysmex Corporation
Kobe, Japan

Georg E. Hoffmann, M.D.
Trillium GmbH
Grafrath, Germany

Stephen Howlett
Beckman Coulter, Inc.
Miami, Florida

Brian Richard Jackson, M.D.
ARUP Laboratories
Salt Lake City, Utah

Paul W. Landesman, Ph.D.
Abbott Laboratories
Abbott Park, Illinois

Michael D. McNeely, M.D.
MDS Metro Clinical Laboratories
Victoria, British Columbia, Canada

Richard A. McPherson, M.D.
Virginia Commonwealth University
Richmond, Virginia

David O'Bryan, Ph.D.
Hibernia Consulting
Kennett Square, Pennsylvania

Paul J. Orsulak, M.D.
VA North Texas Health Care System
Dallas, Texas

Jeff Quint, Ph.D.
Beckman Coulter, Inc.
Brea, California

Richard Seaberg
North Shore University Hospital
Manhasset, New York

Hiroski Sekiya
Olympus America Inc.
Irving, Texas

Russell H. Tomar, M.D.
Cook County Hospital
Chicago, Illinois

Terry Weakley
Cerner Corporation
Kansas City, Missouri

Staff

David E. Sterry, M.T.(ASCP)
Staff Liaison
NCCLS
Wayne, Pennsylvania

Donna M. Wilhelm
Editor
NCCLS
Wayne, Pennsylvania

Melissa A. Lewis
Assistant Editor
NCCLS
Wayne, Pennsylvania

Contents

Abstract.....	i
Committee Membership.....	iii
Foreword.....	ix
1 Scope.....	1
2 Definitions	1
3 Significance and Use.....	2
4 Information Requirements in Clinical Testing.....	3
4.1 General Approach.....	3
4.2 Logical Structure of the Message Level Protocol (See Figure 1.)	4
5 Message Content—General Considerations	5
5.1 Character Codes	5
5.2 Maximum Field Lengths.....	5
5.3 Maximum Record Length.....	6
5.4 Delimiters.....	6
5.5 Data Record Usage Overview.....	8
5.6 Common Field Types.....	9
5.7 Examples of Basic Record Types	11
6 Message Header Record	11
6.1 Record Type ID	11
6.2 Delimiter Definition.....	11
6.3 Message Control ID	11
6.4 Access Password.....	12
6.5 Sender Name or ID	12
6.6 Sender Street Address	12
6.7 Reserved Field	12
6.8 Sender Telephone Number	12
6.9 Characteristics of Sender	12
6.10 Receiver ID	12
6.11 Comment or Special Instructions.....	12
6.12 Processing ID.....	12
6.13 Version Number.....	13
6.14 Date and Time of Message	13
7 Patient Information Record.....	13
7.1 Record Type.....	13
7.2 Sequence Number	13
7.3 Practice-Assigned Patient ID	13
7.4 Laboratory-Assigned Patient ID	13
7.5 Patient ID Number 3	13
7.6 Patient Name.....	13
7.7 Mother's Maiden Name	13
7.8 Birthdate.....	13
7.9 Patient Sex	14

Contents (Continued)

7.10	Patient Race-Ethnic Origin	14
7.11	Patient Address	14
7.12	Reserved Field	14
7.13	Patient Telephone Number	14
7.14	Attending Physician ID.....	14
7.15	Special Field 1	14
7.16	Special Field 2	14
7.17	Patient Height	14
7.18	Patient Weight.....	15
7.19	Patient's Known or Suspected Diagnosis	15
7.20	Patient Active Medications	15
7.21	Patient's Diet.....	15
7.22	Practice Field Number 1	15
7.23	Practice Field Number 2	15
7.24	Admission and Discharge Dates	15
7.25	Admission Status	15
7.26	Location	15
7.27	Nature of Alternative Diagnostic Code and Classifiers	15
7.28	Alternative Diagnostic Code and Classification	16
7.29	Patient Religion.....	16
7.30	Marital Status.....	16
7.31	Isolation Status.....	16
7.32	Language.....	16
7.33	Hospital Service.....	17
7.34	Hospital Institution	17
7.35	Dosage Category.....	17
8	Test Order Record.....	17
8.1	General.....	17
8.2	Multiple Orders.....	17
8.3	General Applications	18
8.4	Field Definitions	18
9	Result Record.....	22
9.1	Record Type ID	22
9.2	Sequence Number	22
9.3	Universal Test ID	22
9.4	Data or Measurement Value	22
9.5	Units.....	23
9.6	Reference Ranges	23
9.7	Result Abnormal Flags	23
9.8	Nature of Abnormality Testing.....	23
9.9	Result Status	23
9.10	Date of Change in Instrument Normative Values or Units	24
9.11	Operator Identification.....	24
9.12	Date/Time Test Started	24
9.13	Date/Time Test Completed.....	24
9.14	Instrument Identification.....	24
10	Comment Record	24
10.1	Record Type ID	24

Contents (Continued)

10.2	Sequence Number	25
10.3	Comment Source.....	25
10.4	Comment Text	25
10.5	Comment Type	25
11	Request Information Record	25
11.1	Record Type ID	25
11.2	Sequence Number	25
11.3	Starting Range ID Number	25
11.4	Ending Range ID Number	26
11.5	Universal Test ID	26
11.6	Nature of Request Time Limits.....	26
11.7	Beginning Request Results Date and Time	26
11.8	Ending Request Results Date and Time.....	26
11.9	Requesting Physician Name	26
11.10	Requesting Physician Telephone Number	27
11.11	User Field Number 1.....	27
11.12	User Field Number 2.....	27
11.13	Request Information Status Codes.....	27
12	Message Terminator Record	27
12.1	Record Type ID	27
12.2	Sequence Number	27
12.3	Termination Code	27
13	Scientific Record.....	28
13.1	Record Type ID	28
13.2	Sequence Number	28
13.3	Analytical Method	28
13.4	Instrumentation	28
13.5	Reagents.....	28
13.6	Units of Measure.....	28
13.7	Quality Control	28
13.8	Specimen Descriptor.....	28
13.9	Reserved Field	28
13.10	Container.....	29
13.11	Specimen ID	29
13.12	Analyte.....	29
13.13	Result	29
13.14	Result Units.....	29
13.15	Collection Date and Time	29
13.16	Result Date and Time	29
13.17	Analytical Preprocessing Steps.....	29
13.18	Patient Diagnosis	29
13.19	Patient Birthdate	29
13.20	Patient Sex	29
13.21	Patient Race	29
14	Manufacturer Information Record	30
	References.....	34

Contents (Continued)

Appendix. Comparison of LIS02 and LIS5	35
Summary of Delegate Comments and Area Committee Responses	36
The Quality System Approach.....	38
Related NCCLS Publications.....	39

Foreword

In 2001, ASTM Committee E31 decided to restructure its operations, with the intent of focusing on standards-development issues such as security, privacy, and the electronic health record. Part of the reorganization plan was to transfer responsibility for E31.13 standards to NCCLS.

Following this transfer, nine standards (formerly ASTM E792; E1029; E1238; E1246; E1381; E1394; E1466; E1639; and E2118) were redesignated as NCCLS standards LIS1 through LIS9. This collection of former ASTM standards provides a wide variety of information relating to clinical laboratory computer systems. Some included documents are of general interest as reference sources; others represent specifications of primary importance to instrument manufacturers. LIS02 is a revision of the former ASTM E1394-97.

The Area Committee on Automation and Informatics has assumed responsibility for maintaining the documents and will revise or update each document in accord with the NCCLS Administrative Procedures. The area committee prioritized LIS2-A as the first standard from this collection to be updated, incorporated into the NCCLS document template, and advanced through the NCCLS consensus process. The area committee will revise other documents in the series in a similar manner.

With the transfer of the former ASTM standards, the Area Committee on Automation and Informatics has expanded its Mission Statement to include laboratory information systems. In the future, the area committee will develop additional standards addressing informatics issues as well as issues related to the integration of patient clinical data.

The revisions in this version of the LIS02 standard are intended to delineate this document from the former ASTM version of this standard. The title and text have been revised throughout to indicate that this standard applies to clinical laboratory instruments. The term computer has been replaced with the term information to better reflect the current terminology (i.e., LIS).

Key Words

Component field, delimiter, field, message, record, repeat field

Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems; Approved Standard—Second Edition

1 Scope

This standard covers the two-way digital transmission of remote requests and results between clinical laboratory instruments and information systems. It is intended to document the common conventions required for the interchange of clinical results and patient data between clinical laboratory instruments and information systems. This standard specifies the message content for transferring information between a clinical laboratory instrument and an information system. It enables any two such systems to establish a logical link for communicating text to send result, request, or demographic information in a standardized and interpretable form. This standard does not necessarily apply to general analytical instruments in an industrial analytical setting, or to a research and development setting.

This standard is intended to apply to the structure of messages exchanged between clinical laboratory instruments and information systems by means of defined communications protocols. Low-level communications protocols and data transfer requirements are beyond the scope of this standard. A separate specification is available detailing a standard for low-level data transfer communications (see NCCLS document LIS1—*Standard Specification for Low-Level Protocol to Transfer Messages Between Clinical Laboratory Instruments and Computer Systems*).

This standard specifies the conventions for structuring the content of the message and for representing the data elements contained within those structures. It is applicable to all text-oriented clinical instrumentation. It has been specifically created to provide common conventions for interfacing computers and instruments in a clinical setting. It would also be applicable to interfacing instruments in clinical practice settings, such as physicians' offices, clinics, and satellite laboratories. The intended users of this standard are developers of clinical laboratory information systems and clinical laboratory managers.

2 Definitions

Battery – A group of tests ordered together, for example, an admitting battery; **NOTES:** a) The term *battery* is used in the document synonymously with the term *profile* or *panel*; b) The test elements within a battery may be characteristic of a single physiologic system, for example, liver function tests, or many different physiologic systems; c) The battery is simply a convention by which a user can order multiple tests by specifying a single name.

Component field – A single data element or data elements which express a finer aggregate or extension of data elements which precede it, for example, parts of a field or repeat field entry; **NOTES:** a) As an example, the patient's name is recorded as last name, first name, and middle initial, each of which is separated by a component delimiter; b) Components cannot contain repeat fields.

Download – Data transmitted from an information system to a clinical instrument.

Field – One specific attribute of a record which may contain aggregates of data elements further refining the basic attribute.

Message – A textual body of information consisting of a header (H) record through a message terminator (L) record.

Record – An aggregate of fields describing one aspect of the complete message.

Repeat field – A single data element which expresses a duplication of the field definition it is repeating;
NOTE: It is used for demographics, requests, orders and the like, where each element of a repeat field is to be treated as having equal priority or standing to associated repeat fields.

Test – A qualitative, semiquantitative, quantitative, or semiquantitative procedure for detecting the presence of, or measuring the quantity of an analyte.

Upload – Data transmitted from a clinical instrument to an information system.

3 Significance and Use

This standard provides for two-way transmission, allowing for data flow in either direction. It provides for sending demographic and test information to or from clinical instruments. This document has sufficient flexibility to permit the addition of fields to existing record types or the creation of new record types to accommodate new test and reporting methodologies.

This document is related to NCCLS document LIS5—*Standard Specification for Transferring Clinical Observations Between Independent Computer Systems*. Both standards use positional convention to define the structure of messages that exchange information about clinical test requests and results. The set of conventions specifies a hierarchical set of records in which the records higher in the hierarchy contain information that is common to all records lower in the hierarchy and thus avoids redundancy in linking data together. The positional convention is simple and direct to implement, requiring only a sequence of strings, each having variable length-delimited fields which are positionally specified.

NCCLS document LIS5—*Standard Specification for Transferring Clinical Observations Between Independent Computer Systems*, in its entirety, is not appropriate for use as a clinical instrument to information system interface. The conventions of NCCLS document LIS5 regarding record types and the organization of data elements within the records have been adhered to as closely as possible to ensure that common data elements defined there and used within instruments are specified as closely as possible. This facilitates the use of this standard consistent with NCCLS document LIS5 in a number of settings. There are three compelling reasons for developing a separate standard which deviates from NCCLS document LIS5.

The scope of NCCLS document LIS5 is specifically targeted to accommodate information transfer between two independent information systems requiring shared patient demographic and test result data. NCCLS document LIS5 contains extensive requirements and limitations, much of which may be of little, if any, use by clinical instrument systems. Further, clinical instruments have test- and instrument-specific requirements outside the scope of NCCLS document LIS5 and, as such, are not available within the existing NCCLS document LIS5.

The structure of NCCLS document LIS5 provides great flexibility in the ordering and reporting of test results and patient demographics. While this is appropriate for use by advanced information systems of equivalent rank, NCCLS document LIS5 clearly falls beyond the technical limitations of many clinical laboratory instruments. This document attempts to identify and simplify all complex data structures and interface procedures and, where practical, restrict multiple procedural options to single procedures appropriate for the clinical instrument setting. Further, this document has attempted to assign a master/slave hierarchy where conflicts may occur, assigning appropriate responsibility for data handling or reporting operations to the party (clinical instrument or information system) better able to process a particular task. For example, in all cases involving the ordering or reporting of tests, the instrument manufacturer is solely responsible for assigning the test and result ID numbers (see Section 5.6.1). These reductions in flexibility directly result in increased structure and clarity, which is deemed more appropriate for ensuring successful interface implementation within the clinical instrument setting.

NCCLS document LIS5 was developed independent of data protocol and transfer considerations. NCCLS document LIS5 uses maximum field and record lengths. Combined with its record level checksum and error recovery facilities, NCCLS document LIS5 may be implemented without a data protocol layer. By contrast, this message-content specification has been developed in cooperative effort with a correlative low-level data transfer and protocol specification. While each specification (message-content and low-level protocol) is designed to be independently implemented and maintained, the message-content specification presumes that a protocol layer exists that will handle record blocking/deblocking, error detection and recovery, and other associated data transport tasks. As such, all protocol level operations and limitations existing in NCCLS document LIS5 are not applicable, and therefore not included in this document.

4 Information Requirements in Clinical Testing

4.1 General Approach

Messages may contain one or more requests/results for one or more patients. Tests may be requested as groups of many individual tests. These groups are referred to as “batteries.” Examples of batteries are tests produced on a multichannel analyzer, such as a basic metabolic profile (BMP), physiological groupings of tests (such as liver function tests) and minimal inhibitory concentration tests (MICs) in microbiology testing. The fact that a series of tests is contained in a battery does not imply that they are all performed on the same analytic instrument.

Messages consist of a hierarchy of records of various types. Records at level zero contain information pertaining to the sender identification and completion of transmission. Records at level one of the hierarchy contain information about individual patients. Records at level two contain information about test order requests and specimens. Records at level three contain information about test results.

Comment records may be inserted at any level in the hierarchy. A comment record always relates to the immediately preceding patient, order, result, scientific, or manufacturer information record. Therefore, if a comment record were to follow a patient record (level one), then that comment record would be treated as a level two record. A comment record may not follow the message terminator record.

Manufacturer information records may be inserted at any level in the hierarchy. This record type always relates to the immediately preceding patient, order, result, scientific, or comment record. Therefore, if a manufacturer information record were to follow a patient record (level one), then the record would be treated as a level two record. This record may not follow the message terminator record.

Additional record types are the request-information record and the terminator record. The request-information record provides for the request of demographics or test results to or from the clinical instrument for specified patients, specimens, tests, and dates, and the like. The message terminator record must be the very last record of the message.

The smallest element of information in any record is the field, containing a single item of information, such as a date, a patient name, or a numeric test result.

The test order record contains information about ordering a single test, test battery, or a series of tests or batteries, as discussed in Sections 5.5.3 and 8.

Most of the record types are related to each other in a definite hierarchy. At level zero is the message header and message terminator. At level one is the patient record, the request-information record, and the scientific record. At level two is the test order record. At level three is the result record. The comment and manufacturer information records do not have an assigned level.

A sequence of patient records, order records, or result records at one level is terminated by the appearance of a record type of the same or higher level. Thus, a sequence of results for one battery of tests is terminated by the next test order, patient, manufacturer information, request information, or message terminator record.

An order record may never appear without a preceding patient record, and a result record may never appear without a preceding order record.

When an order is transmitted, it must be preceded by a patient record. All orders that follow apply to the patient in the preceding patient record. When a result is transmitted, it must be preceded by an order record and a patient record to maintain the prescribed hierarchy.

Each instrument manufacturer adhering to this standard may decide which fields are applicable for their particular application with the exception of those fields necessary to identify the record type or parse individual fields. Thus, the need to send the hierarchy of records need not generate large messages.

4.2 Logical Structure of the Message Level Protocol (See Figure 1.)

4.2.1 Logical Information Storage Requirements

In order to determine buffering requirements, both transmitter and receiver must use common rules for storing transmitted data in order to ensure proper error logging and error recovery procedures (see Section 4.2.2). Since data content is structured in a hierarchical fashion, any decremental change in the hierarchical level shall trigger storage of all data transmitted prior to said level change. This rule may be considered as the minimal implementation. Data may be saved at more frequent intervals at the receiver's option (see Figure 2).

4.2.2 Logical Transmission Error Recovery Requirements

Transmission line failure, determined at the transmission protocol level, requires a mechanism for restarting the incomplete message. If a transmission failure occurs, transmission shall restart at the last logical record not presumed saved as outlined in Section 4.2.1. Procedures for determining time before retransmission or maximum number of retransmissions are not within the scope of this document. In order to fulfill hierarchical record level requirements, all logical records necessary to reach the restart record point must be repeated prior to transmitting the record where line failure originally occurred. Using the transmission example as given in Section 4.2.1 and as outlined in Figure 2, the following record recovery examples would be valid.

Line Failure Occurs At:	Requires Retransmission Of:
A	A
B	A, B
C	A, B, C
D	A, B, C, D
E	A, B, C, D, E
F	A, B, E, F
G	A, B, E, F, G
H	A, G, H
I	A, G, H, I
J	A, G, H, I, J
K	A, G, H, I, J, K
L	A, G, H, I, J, K, L
M	A, G, H, L, M
N	A, G, M, N
O	A, N, O
P	A, N, O, P
Q	A, N, O, P, Q

5 Message Content—General Considerations

5.1 Character Codes

All data shall be represented as eight-bit, single-byte, coded graphic character values as defined in ISO 8859-1:1987.¹ The eight-bit values, within the range from 0 to 127 of ISO 8859-1:1987¹ correspond to the ASCII standard character set. Values from 0 to 31 are disallowed with the exception of 7 (BEL), 9 (Horizontal tab), 11 (Vertical tab), and 13 (CR), where 13 is reserved as a record terminator. Values from 32 to 126 and from 128 to 254 are allowed. Values 127 and 255 are also not allowed. It is the responsibility of the instrument vendor and information system vendor to understand the representation of any extended or alternate character set being used. As an example, the numeric value 13.5 would be sent as four-byte value characters 13.5 or Latin-1(49), Latin-1(51), Latin-1(46), and Latin-1(53).

Allowed Characters: 7, 9, 11, 12, 13, 32 to 126, 128 to 254

Disallowed Characters: 0 to 6, 8, 10, 14 to 31, 127, 255

Within text data fields, only the Latin-1 characters 32 to 126 and the undefined characters 128 to 254 are permitted as usable characters (excluding those used as delimiter characters in a particular transmission). Furthermore, all characters used as delimiters in a particular transmission are excluded from the permitted range. The sender is responsible for screening all text data fields to ensure that the text does not contain those delimiters. Unless otherwise stated, contents of data fields shall be case sensitive.

5.2 Maximum Field Lengths

This specification assumes that all fields are variable in length. No storage is allocated (except for the delimiter) for a null field. When, for example, ten characters of data are entered within a field, only ten characters will be used. This specification does not define a maximum length for any field or record and relies upon the receiver's buffering capabilities, and the logical layer's transport facilities, to parse information into workable lengths for transmission and processing purposes. It is the responsibility of the instrument vendor and information system vendor to agree on any arbitrary field or record truncation that may need to be imposed. It is recommended that the instrument vendor provide documentation disclosing any field or record limits that will be mandated by the clinical instrument.

5.3 Maximum Record Length

There is no imposed maximum record length.

5.4 Delimiters

Alphanumeric characters should not be used as delimiters, because they are likely to appear within field content. Moreover, some alphabetic characters have special uses as follows:

H, P, O, R, C, Q, S, L, M	record type IDs
. Latin-1(46)	decimal point (period)
, Latin-1(44)	comma
S, P, R, C	priority codes
L, H, <, >, N, U, D, B, W	result codes
C, P, F, X, I, O	result status

For the purpose of providing examples, the following delimiters are used in this specification:

Field delimiter = vertical bar (|) Latin-1(124)
 Repeat delimiter = backslash (\) Latin-1(96)
 Component delimiter = caret (^) Latin-1(94)
 Escape delimiter = ampersand (&) Latin-1(38)

5.4.1 Record Delimiter

Carriage return (13) shall be the delimiter for the end of any of the defined record types.

5.4.2 Field Delimiter

A single allowable character as defined in Section 5.1, excluding Latin-1(13) (carriage return), shall separate adjacent fields. The field delimiter is variable and defined in the message header. The same delimiter must be used in all records following a header and preceding a message terminator record.

5.4.3 Repeat Delimiter

A repeat delimiter is a single allowable character as defined in Section 5.1, excluding Latin-1(13) and the value for the field delimiter defined in Section 5.4.2. The repeat delimiter must be defined in the message header and is used to separate variable numbers of descriptors for fields containing parts of equal members of the same set.

5.4.4 Component Delimiter

A component delimiter is a single allowable character as defined in Section 5.1, excluding Latin-1(13) and the field and repeat delimiter values. The component delimiter is used to separate data elements of fields of a hierarchical or qualifier nature. For example, the street, city, state, zip, etc. of an address field would be separated by component delimiters.

5.4.5 Escape Delimiter

An escape delimiter is a single allowable character, as defined in Section 5.1, excluding Latin-1(13) and the field, repeat, and component delimiter values. The escape delimiter is used within text fields to signify special case operations. Applications of the escape delimiter are optional and may be used or ignored at

the discretion of either the transmitter or the receiver. However, all applications are required to accept the escape delimiter and use it to correctly parse fields within the record.

5.4.5.1 Use of Escape Delimiter

The escape delimiter may be used to signal certain special characteristics of portions of a text field (e.g., imbedded delimiters, line feed, carriage return). An escape sequence consists of the escape delimiter character followed by a single escape code ID (listed below), followed by zero or more data characters followed by another (closing) occurrence of the escape delimiter character. No escape sequence may contain a nested escape sequence. The following escape sequences are predefined.

&H&	start highlighting text
&N&	normal text (end highlighting)
&F&	imbedded field delimiter character
&S&	imbedded component field delimiter character
&R&	imbedded repeat field delimiter character
&E&	imbedded escape delimiter character
&Xhhhh&	hexadecimal data

NOTE: Any number of hexadecimal digits (0 to 9, A to F) may follow (that is, &XA& could equal line feed).

&Zcccc&	Local (manufacturer-defined) escape sequence
---------	--

NOTE: Any number of legal characters may follow.

5.4.6 Specification of Delimiters

The actual delimiters to be employed in a given transmission shall be specified in the header message. It is the responsibility of the sender to avoid the inclusion of any delimiter characters within the field contents. The receiving information system will determine what characters to use by reading the specifications of the header it receives. See Section 5.4 for examples of delimiters used for this document.

5.4.7 Delimiters for Null Values

Fields shall be identified by their position, obtained by counting field delimiters from the front of the record. This position-sensitive identification procedure requires that when the contents of the field are null, its corresponding field delimiter must be included in the record to ensure that the *i*'th field can be found by counting (*i*–1) delimiters. Delimiters are not included for trailing null fields; that is, if the tenth field was the last field containing data, the record could terminate after the tenth field, and therefore would contain only nine delimiters.

5.4.8 Fields of No Concern to the Receiving System

Transmitted records may include more fields than are required by a receiving system. When processing a message, the receiving system may ignore any field it does not require. Fields must always be transmitted, however, in the positional order specified.

5.4.9 Fields with Null Values

A system may transmit a null value for a field because: 1) it does not know the value; 2) it knows the value is irrelevant to the receiving system; or 3) the value has not changed since the last transmission, or any combination thereof. To exemplify the third case, a laboratory within a tightly linked hospital

network may never transmit the patient's birthdate, sex, or race in the patient record when transmitting the order and result records to the requesting system, because it knows that the hospital registry system always broadcasts new or changed patient data to the receiving system.

Because the sending system can use null values to indicate no change, a null value does not overwrite existing data in the receiving system. In rare circumstances, for example, if a system erroneously sent a patient's birthdate when the birthdate was actually unknown, the receiving system should replace its existing value for a field with a null value.

A field containing only a pair of double quotes (Latin-1(34)) should be treated as an instruction to the receiver that the existing contents pertaining to that field definition should be deleted.

The use of null fields to default to previous values is discouraged, especially in the case of specimen numbers, where it can cause significant ambiguity. This section is retained for compatibility with previous versions of this standard. The vendor must document the behavior of the handling of null fields by the interface.

5.5 Data Record Usage Overview

Data shall be exchanged in records of different types. Each record is introduced by field (number one) identifying the record type, and terminated by a carriage return. The following record types are defined.

NOTE: The record type ID field shall be case insensitive.

5.5.1 Message Header Record (H)

This record shall contain information about the sender and the receiver, that is, it shall identify the instrument(s) and the information systems whose records are being exchanged. It also defines the field, repeat field, and component field delimiter characters.

5.5.2 Patient Identifying Record (P)

This record type contains information about an individual patient.

5.5.3 Test Order Record (O)

When sent from the information system to the instrument, this record shall represent a test order and may be followed by one or more result records which would contain information pertinent to the test being ordered. When sent by the instrument to the information system, it shall provide information about the specimen/test request, and may be followed by result records (at least one record for each test within the ordered batteries).

5.5.4 Result Record (R)

Each result record shall contain the results of a single analytic determination.

5.5.5 Comment Record (C)

Comment records shall apply to any other record except the message trailer record. They may be free standing messages sent to or from the instrument, unrelated to a particular patient or test procedure.

5.5.6 Request Information Record (Q)

This record shall be used to request information for new tests, tests previously ordered, and possibly tests previously reported. A single request information record may request demographic information, or results for an individual test, multiple tests, or all tests for a single date, a series of dates, or a range of dates, or both, and for an individual patient, group of patients, individual specimens, groups of specimens, etc.

5.5.7 Scientific Record (S)

This record shall be used to exchange results between clinical sites for the purposes of proficiency testing or method development.

5.5.8 Manufacturer Information Record (M)

This record, which is similar to the comment record, may be used to send complex structures where use of the existing record types would not be appropriate. The fields within this record type are defined by the manufacturer.

5.6 Common Field Types

5.6.1 Universal Test ID

This field is defined as a four-part field with provisions to further define the test identification via use of component fields. The test ID field is used to identify a test or battery name. The four parts which are defined below are the universal test identifier, the test name, the test identifier type, and the manufacturer-defined test code. All test ID parts must be separated by a component delimiter and are position dependent. As an example, additional information which may be included in this field type are instrument ID, organism ID (for sensitivity tests), well number, cup number, location number, tray number, bar-code number, etc. It is the responsibility of the instrument manufacturer to define the data content of the test ID field. When the test ID is used in the result record, there must be sufficient information within the test ID field to determine the relationship of the test result to the test, battery, or batteries ordered.

5.6.1.1 Universal Test ID (Part 1)

This is the first component of the test ID field. This field is currently unused but reserved for the application of a universal test identifier code (LOINC Codes), should one system become available for use at a future time.

5.6.1.2 Universal Test ID Name (Part 2)

This would be the test or battery name associated with the universal test ID code described in Section 5.6.1.1.

5.6.1.3 Universal Test ID Type (Part 3)

In the case where multiple national or international coding schemes exist, this field may be used to determine what coding scheme is employed in the test ID and test ID name fields.

5.6.1.4 Manufacturer's or Local Code (Part 4)

This is the code defined by the manufacturer. This code may be a number, characters, or a multiple test designator based on manufacturer-defined delimiters (that is, AK.23.34-B). Extensions or qualifiers to this code may be followed by subsequent component fields which must be defined and documented by the

manufacturer. For example, this code may represent a three-part identifier such as—Dilution ^ Diluent ^ Description.

5.6.2 Dates and Times

In all cases, dates shall be recorded in the YYYYMMDD format as required by ANSI X3.30.² December 1, 1989 would be represented as 19891201. When times are transmitted, they shall be represented as HHMMSS, and shall be linked to dates as specified by ANSI X3.43.³ Date and time together shall be specified as up to a 14-character string: YYYYMMDDHHMMSS.

5.6.2.1 Time Zone

The time zone may be optionally appended to the date/time field in the format + HHMM or – HHMM as appropriate. The default time zone is that of the sender.

5.6.3 Telephone Numbers

Phone numbers shall be recorded as free text, which may contain extensions such as area code, country code, beeper number, and hours to call.

5.6.3.1 Multiple Phone Numbers

When multiple telephone numbers apply, they may be included in one field and separated from each other by repeat delimiters. The first such entry is considered the primary or the daytime number.

5.6.4 Fixed Measurements and Units

When a field contains a specific observation, for example, patient's weight, patient's height, or collection volume, the default units of measurement for that observation are specified in the field definition. When the observation is measured in the default units, the units need not be transmitted. If the measure is recorded in units different from the default, for example, if the weight is measured in pounds rather than kilograms, the measurement units must be transmitted. In this case, the units are transmitted in the same field as the measurement. The units follow the measure and are separated from it by a component delimiter, for example, 100 ^ lb. Units should be expressed in ISO standard abbreviations in accordance with ISO 2955.⁴

5.6.5 Addresses

An address occupies a single field in a record. The address may be comprised of five components (street address, city, state, zip, or postal code, and country code) separated by component delimiters so that the receiving party can break them into separate fields as needed. An example would be 52 Hilton Street #B42 ^ Chicago ^ IL ^ 60305 ^ USA. The country need only be transmitted when it cannot be assumed from the context. The components of this field are position dependent.

5.6.6 Provider and User IDs

Physicians' and other caregivers' codes may be transmitted as internal code numbers, as full names, or both, as mutually agreed upon between the sender and the receiver. When both the name and ID number are sent, ID numbers should come first and be separated from the name by a component delimiter. Each component of the name is also separated by a component delimiter. The order of the components of the name shall be: 1) last name; 2) first name; 3) middle initial or name; 4) suffix (e.g., Jr., Sr.); and 5) title (e.g., Dr., Mr.). Thus, if Dr. John G. Jones, Jr. had an identifier of 401-0, his number and name would be

transmitted as 401-0 ^ JONES ^ JOHN ^ G ^ JR ^ DR. When necessary, more than one ID may be sent within one field. Multiple IDs in one field are separated by repeat delimiters.

5.6.7 Record Sequence Number

This is a *required* field used in record types that may occur multiple times within a single message. The number used defines the *i*'th occurrence of the associated record type at a particular hierarchical level and is reset to one whenever a record of a greater hierarchical significance (lower number) is transmitted or if the same record is used at a different hierarchical level (e.g., comment records).

5.7 Examples of Basic Record Types

The following examples are given for a set of transmitted results for a given patient. These will show how the employment of the conventions defined lead to a valid message. In these examples, the first two fields of each line (record) of the message body contain the record type and the integer record sequence number (excepting the header record). Carriage return is indicated by <CR>. To simplify the example, all the components of each record have not been included. Ellipses (...) are used to indicate fields that are left out and comments are enclosed in square brackets. Record hierarchical levels are shown by indentation.

NOTE: The user may wish to study the record definitions outlined in Section 6 before reviewing the samples shown in Figures 3 to 11. Trailing fields, unused, may or may not have field delimiters transmitted. Both cases should be handled by the receiving parser.

6 Message Header Record

The header shall contain identifiers of both the sender and the receiver. The message header is a level zero record and must be followed at some point by a message terminator record before ending the session or transmitting another header record. This record type must always be the first record in a transmission.

6.1 Record Type ID

The character H identifies the record as a message header record.

6.2 Delimiter Definition

The five Latin-1 characters that immediately follow the H (the header ID) define the delimiters to be used throughout the subsequent records of the message. The second character in the header record is the field delimiter, the third character is the repeat delimiter, the fourth character is the component delimiter, and the fifth is the escape character. A field delimiter follows these characters to separate them from subsequent fields. Another way to view this is that the first field contains H and the second field contains the repeat, component, and escape delimiters. Using the example delimiters, the first six characters in the header record would appear as follows: H | ^ & |.

6.3 Message Control ID

This is a unique number or other ID that uniquely identifies the transmission for use in network systems that have defined acknowledgment protocols that are outside of the scope of this document. Note that this is the third field.

6.4 Access Password

This is a level security/access password as mutually agreed upon by the sender and receiver. If this security check fails, the transmission will be aborted, and the sender will be notified of an access violation.

6.5 Sender Name or ID

The purpose of this field is to define the manufacturer/instrument(s) specific to this line. Using repeat and/or component delimiters, this field may reflect software or firmware revisions, multiple instruments available on the line, etc.

6.6 Sender Street Address

This text value shall contain the street address of the sender as specified in Section 5.6.5.

6.7 Reserved Field

This field is currently unused but reserved for future use.

6.8 Sender Telephone Number

This field identifies a telephone number for voice communication with the sender as specified in Section 5.6.3.

6.9 Characteristics of Sender

This field contains any characteristics of the sender such as, parity, checksums, optional protocols, etc. necessary for establishing a communication link with the sender.

6.10 Receiver ID

This text value includes the name or other ID of the receiver. Its purpose is verification that the transmission is indeed for the receiver.

6.11 Comment or Special Instructions

This text field shall contain any comments or special instructions relating to the subsequent records to be transmitted.

6.12 Processing ID

The processing ID indicates how this message is to be processed:

- P Production: Treat message as active message to be completed according to standard processing.
- T Training: Message is initiated by a trainer and should not have an effect on the system.
- D Debugging: Message is initiated for the purpose of a debugging program.
- Q Quality Control: Message is initiated for the purpose of transmitting quality control/quality assurance or regulatory data.

6.13 Version Number

This value identifies the version level of the specification. This value is currently LIS02-A2.

6.14 Date and Time of Message

This field contains the date and time that the message was generated using the format specified in Section 5.6.2.

7 Patient Information Record

Each line of the patient record shall begin with a record type and end with a carriage return.

7.1 Record Type

The character P identifies the record as a patient record.

7.2 Sequence Number

For the first patient transmitted, 1 shall be entered, for the second, 2, ... until the last as defined in Section 5.6.7.

7.3 Practice-Assigned Patient ID

This identifier shall be the unique ID assigned and used by the practice to identify the patient and his/her results upon return of the results of testing.

7.4 Laboratory-Assigned Patient ID

This identifier shall be the unique processing number assigned to the patient by the laboratory.

7.5 Patient ID Number 3

This field shall be optionally used for additional, universal, or manufacturer-defined identifiers (such as the social security account no.), as arranged between the transmitter and the receiver. Please note that individuals are not required to provide social security numbers.

7.6 Patient Name

The patient's name shall be presented in the following format: last name, first name, middle name or initial, suffix, and title, and each of these components shall be separated by a component delimiter as described in Section 5.6.6.

7.7 Mother's Maiden Name

The optional mother's maiden name may be required to distinguish between patients with the same birthdate and last name when registry files are very large. This name shall be presented as the mother's maiden surname, for example, Thompson.

7.8 Birthdate

The birthdate shall be presented in the standard format specified in Section 5.6.2.

7.9 Patient Sex

This field shall be represented by M, F, or U.

7.10 Patient Race-Ethnic Origin

The following examples may be used:

W	white
B	black
O	Asian/Pacific Islander
NA	Native American/Alaskan Native
H	Hispanic

Full text names of other ethnic groups may also be entered. Note that multiple answers are permissible, separated by a component delimiter.

7.11 Patient Address

This text value shall record the street address of the patient's mailing address as defined in Section 5.6.5.

7.12 Reserved Field

This field is reserved for future expansion.

7.13 Patient Telephone Number

The patient's telephone number is formatted as defined in Section 5.6.3.

7.14 Attending Physician ID

This field shall identify the physician(s) caring for the patient as either names or codes, as agreed upon between the sender and the receiver. Identifiers or names, or both, should be separated by component delimiters as specified in Section 5.6.6. Multiple physician names (for example, ordering physician, attending physician, referring physician) shall be separated by repeat delimiters.

7.15 Special Field 1

This is an optional text field for vendor use (each laboratory can use this differently).

7.16 Special Field 2

This is an optional text field for vendor use.

7.17 Patient Height

This is an optional numeric field containing the patient's height. The default units are centimeters. If measured in terms of another unit, the units should also be transmitted as specified in Section 5.6.4.

7.18 Patient Weight

This is an optional numeric field containing the patient's weight. The default units are kilograms. If measured in terms of another unit, for example, pounds, the unit name shall also be transmitted as specified in Section 5.6.4. Height and weight information is not currently required by all laboratories but is of value in estimating normative values based upon body surface area.

7.19 Patient's Known or Suspected Diagnosis

This value should be entered either as an ICD-9 code or as free text. If multiple diagnoses are recorded, they shall be separated by repeat delimiters.

7.20 Patient Active Medications

This field is used for patient active medications or those suspected, in overdose situations. The generic name shall be used. This field is of use in interpretation of clinical results.

7.21 Patient's Diet

This optional field in free text should be used to indicate such conditions that affect results of testing, such as 16-hour fast (for triglycerides) and no red meat (for hemocult testing).

7.22 Practice Field Number 1

This is a text field for use by the practice; the optional transmitted text will be returned with the results.

7.23 Practice Field Number 2

This is the same as described in Section 7.22.

7.24 Admission and Discharge Dates

These values shall be represented as specified in Section 5.1. The discharge date, when included, follows the admission date and is separated from it by a repeat delimiter.

7.25 Admission Status

This value shall be represented by the following minimal list or by extensions agreed upon between the sender and receiver: OP (outpatient), PA (preadmit), IP (inpatient), ER (emergency room).

7.26 Location

This text value shall reflect the general clinic location or nursing unit, or ward or bed (or both) of the patient in terms agreed upon by the sender and the receiver.

7.27 Nature of Alternative Diagnostic Code and Classifiers

This field relates to Section 7.28. It identifies the class of code or classifiers that are transmitted (e.g., DRGs, or in the future, AVGs [ambulatory visitation groups]).

7.28 Alternative Diagnostic Code and Classification

Alternative diagnostic codes and classifications (e.g., DRG codes) can be included in this field. The nature of the diagnostic code is identified in Section 7.27. If multiple codes are included, they should be separated by repeat delimiters. Individual codes can be followed by optional test descriptors (when the latter are present) and must be separated by component delimiters.

7.29 Patient Religion

When needed, this value shall include the patient's religion. Codes or names may be sent as agreed upon between the sender and the receiver. Full names of religions may also be sent as required. A list of sample religious codes follows:

P	Protestant
C	Catholic
M	Church of the Latter Day Saints (Mormon)
J	Jewish
L	Lutheran
H	Hindu

7.30 Marital Status

When required, this value shall indicate the marital status of the patient as follows:

M	married
S	single
D	divorced
W	widowed
A	separated

7.31 Isolation Status

Isolation codes indicate precautions that must be applied to protect the patient or staff against infection. The following are suggested codes for common precaution. Multiple precautions can be listed when separated by repeat delimiters. Full text precautions may also be sent.

ARP	antibiotic resistance precautions
BP	blood and needle precautions
ENP	enteric precautions
NP	precautions for neutropenic patient
PWP	precautions for pregnant women
RI	respiratory isolation
SE	secretion/excretion precautions
SI	strict isolation
WSP	wound and skin precautions

7.32 Language

The value of this field indicates the patient's primary language. This may be needed when the patient is not fluent in the local language.

7.33 Hospital Service

This value indicates the hospital service currently assigned to the patient. Both code and text may be sent when separated by a component delimiter as in Section 5.6.6.

7.34 Hospital Institution

This value indicates the hospital institution currently assigned to the patient. Both code and text may be sent when separated by a component delimiter as in Section 5.6.6.

7.35 Dosage Category

This value indicates the patient dosage group. For example, A–ADULT, P1–PEDIATRIC (one to six months), P2–PEDIATRIC (six months to three years), etc. Subcomponents of this field may be used to define dosage subgroups.

8 Test Order Record

8.1 General

The test order record defines the attributes of a particular request for a clinical instrument's services and contains all specimen information. An order record will be generated by the information system to request a given test, battery, or set of tests. The information in an order record will usually apply to a single specimen. However, there is not necessarily a one-to-one relationship between specimen and tests ordered. Different test batteries will usually be ordered within different order records even when they can be performed on a single specimen. In this case, the specimen information is duplicated in each of the order records that employ that specimen.

8.2 Multiple Orders

More than one test or test battery may be ordered on a single order record by using repeat delimiters between the individual tests ordered in that record. However, in such cases, all other attributes stored within the order record must be the same for all the tests ordered within that record. Thus, if one wishes to order one test as a STAT or immediate test and another as a routine test, two separate order records would be required. In the case that a test battery requires more than one specimen, such as is true for creatinine clearances, information about each of the test specimens may be included in the single order record identifying multiple specimens using the repeat delimiter within the specimen ID field.

Although multiple tests or test batteries can be on a single order record, when reporting the results, the instrument shall produce a separate order record for each unique battery, copying the appropriate specimen information from the original order record into each of the new order records.

In the event that a test battery cannot be performed, for example, because of hemolysis, the order record will be returned to the information system with the report type indicator X to indicate that it was not done. In this case, no result records will be transmitted.

When test analyses are successfully performed, the message returned to the information system will include the order record followed by result records for each separate observation requested by that order. The number of such result records will depend upon the number of individual measurements performed in the analysis. Four test result records would follow the order record for an electrolytes test. Twelve result records will follow the order record for a BMP.

Test batteries that require multiple specimens for their performance would similarly be followed by a series of result records corresponding to the number of individual measurements obtained. The manufacturer must ensure that the test ID field within each result record contains sufficient information to relate the individual test measurements to the specific tests, batteries, and specimens ordered.

Microbiological culture results are different. A new order record should be created for each panel of antimicrobial sensitivities, although multiple batteries/panels may be ordered on a single order record if desired. The series of antimicrobial sensitivities for any single sensitivity analysis will be reported as separate result records, one for each result element or combination of elements (antimicrobial, MIC, interpretation, etc.). Thus, the antimicrobial sensitivity appears logically very much like an extended BMP result with separate result records for each separate result from each antibiotic tested. Once again, the test ID field within the result records must contain sufficient information to relate the individual test measurements with the appropriate antibiotic test and battery ordered.

8.3 General Applications

The order record may be used in four different circumstances:

- (1) It is sent by the information system to request a particular set of instrument tests.
- (2) It is transmitted back to the information system as part of the results. If the ordered instrument analyses can be completed, the instrument sends back the order record along with the result records according to the hierarchy described in this document. If results cannot be produced, for example, because the specimen is hemolyzed, the laboratory transmits the order with an appropriate report type (see Section 8.4.26) to indicate this problem, but no result records are transmitted.
- (3) The order record is transmitted back to the information system in response to a request information query. In this case, it has the same form as in paragraph (2) above.
- (4) The instrument is requesting demographic or tests-ordered information from the information system.

8.4 Field Definitions

The order record is comprised of the following fields:

8.4.1 Record Type ID

The character assigned to the order record shall be O.

8.4.2 Sequence Number

This field shall be represented as described in Section 5.6.7.

8.4.3 Specimen ID

This text field shall represent a unique identifier for the specimen assigned by the information system and returned by the instrument. If the specimen has multiple components further identifying cultures derived from it, these component identifiers will follow the specimen ID and be separated by component delimiters. For example, the specimen ID may contain the specimen number followed by the isolate number, well or cup number (for example, 10435A ^ 01 ^ 64).

8.4.4 Instrument Specimen ID

This text field shall represent a unique identifier assigned by the instrument, if different from the information system identifier, and returned with results for use in referring to any results.

8.4.5 Universal Test ID

This field shall use universal test ID as described in Section 5.6.1.

8.4.6 Priority

Test priority codes are as follows:

S	stat
A	as soon as possible
R	routine
C	callback
P	preoperative

If more than one priority code applies, they must be separated by repeat delimiters.

8.4.7 Requested/Ordered Date and Time

The contents of this field shall be represented as specified in Section 5.6.2 and will denote the date and time the test order should be considered ordered. Usually this will be the date and time the order was recorded. This is the date and time against which the priorities should be considered. If the ordering service wants the test performed at a specified time in the future, for example, a test to be drawn two days in the future at 8 p.m., the future date and time should be recorded here. Note that the message header data and the future date and time should be recorded here. Further, note that the message header record date and time (see Section 6.14) indicates the time the order was transmitted to or from the instrument.

8.4.8 Specimen Collection Date and Time

This field shall represent the actual time the specimen was collected or obtained.

8.4.9 Collection End Time

This field shall contain the end date and time of a timed specimen collection, such as 24-hour urine collection. The value shall be specified according to Section 5.6.2.

8.4.10 Collection Volume

This value shall represent the total volume of specimens such as urine or other bulk collections when only aliquot is sent to the instrument. The default unit of measure is milliliters. When units are explicitly represented, they should be separated from the numeric value by a component delimiter, for example, 300 ^ g. Units should follow the conventions given in Section 5.6.4.

8.4.11 Collector ID

This field shall identify the person and facility which collected the specimen. If there are questions relating to circumstances surrounding the specimen collection, this person will be contacted.

8.4.12 Action Code

This field shall indicate the action to be taken with respect to the specimens that accompany or precede this request. The following codes shall be used:

C	cancel request for the battery or tests named
A	add the requested tests or batteries to the existing specimen with the patient and specimen identifiers and date/time given in this record
N	new requests accompanying a new specimen
P	pending specimen
L	reserved
X	specimen or test already in process
Q	treat specimen as a Q/C test specimen

8.4.13 Danger Code

This field representing either a test or a code shall indicate any special hazard associated with the specimen, for example, a hepatitis patient, suspected anthrax.

8.4.14 Relevant Clinical Information

Additional information about the specimen would be provided here and used to report information such as amount of inspired O₂ for blood gases, point in menstrual cycle for cervical pap tests, or other conditions that influence test interpretations.

8.4.15 Date/Time Specimen Received

This optional field shall contain the actual log-in time recorded in the laboratory. The convention specified in Section 5.6.2 shall be used.

8.4.16 Specimen Descriptor

This field may contain two separate elements—specimen type and specimen source—as defined in Sections 8.4.16.1 and 8.4.16.2. The components must be separated by component delimiters.

8.4.16.1 Specimen Type

Samples of specimen culture types or sources would be blood, urine, serum, hair, wound, biopsy, sputum, etc.

8.4.16.2 Specimen Source

This is always the second component of the specimen descriptor field and is used specifically to determine the specimen source body site (e.g., left arm, left hand, right lung).

8.4.17 Ordering Physician

This field shall contain the name of the ordering physician in the format outlined in Section 5.6.6.

8.4.18 Physician's Telephone Number

This field shall contain the telephone number of the requesting physician and will be used in responding to callback orders and for critically abnormal results. Use the format given in Section 5.6.3.

8.4.19 User Field Number 1

Text sent by the requestor should be returned by the sender along with the response.

8.4.20 User Field Number 2

This field is similar to that described in Section 8.4.19.

8.4.21 Laboratory Field Number 1

This is an optional field definable for any use by the laboratory.

8.4.22 Laboratory Field Number 2

This field is similar to that described in Section 8.4.21.

8.4.23 Date/Time Results Reported or Last Modified

This field is used to indicate the date and time the results for the order are composed into a report, or into this message or when a status as defined in Sections 8.4.26 or 9.9 is entered or changed. When the information system queries the instrument for untransmitted results, the information in this field may be used to control processing on the communications link. Usually, the ordering service would only want those results for which the reporting date and time is greater than the date and time the inquiring system last received results. Dates and times should be recorded as specified in Section 5.6.2.

8.4.24 Instrument Charge to Information System

This field contains the billing charge or accounting reference by this instrument for tests performed.

8.4.25 Instrument Section ID

This identifier may denote the section of the instrument where the test was performed. In the case where multiple instruments are on a single line or a test was moved from one instrument to another, this field will show which instrument or section of an instrument performed the test.

8.4.26 Report Types

The following codes shall be used:

O	order record; user asking that analysis be performed
C	correction of previously transmitted results
P	preliminary results
F	final results
X	order cannot be done, order cancelled
I	in instrument pending
Y	no order on record for this test (in response to query)
Z	no record of this patient (in response to query)
Q	response to query (this record is a response to a request-information query)

8.4.27 Reserved Field

This field is unused but reserved for future expansion.

8.4.28 Location of Specimen Collection

This field defines the location of specimen collection if different from the patient location.

8.4.29 Nosocomial Infection Flag

This field is used for epidemiological reporting purposes and will show whether the organism identified is the result of a nosocomial (hospital-acquired) infection.

8.4.30 Specimen Service

In cases where an individual service may apply to the specimen collected, and the service is different from the patient record service, this field may be used to define the specific service responsible for such collection.

8.4.31 Specimen Institution

In cases where the specimen may have been collected in an institution, and the institution is different from the patient record institution, this field may be used to record the institution of specimen collection.

9 Result Record

The result record shall include the following fields:

9.1 Record Type ID

The record type ID is coded as R.

9.2 Sequence Number

The sequence number shall be assigned as described in Section 5.6.7.

9.3 Universal Test ID

This field shall use the universal test IDs described in Section 5.6.1.

9.4 Data or Measurement Value

Whether numeric, text, or coded values, the data shall be recorded in ASCII text notation. If the data result contains qualifying elements of equal stature, these should be separated by component delimiters. This applies strictly to results of identical nature (that is, this field may not contain implied subvalues). Use of components within this field should be avoided whenever possible.

Multiple results or values, observed, calculated, or implied, for a single test order (i.e., MIC or interpretation codes from a single antibiotic sensitivity test) must be reported in separate result records with each result definition defined uniquely by the test ID field as given in Section 9.3. Correspondingly, the test ID field (see Section 9.3) must be sufficiently descriptive to determine the placement of the data value with reference to the original test order record and to other result records associated with said test order record.

9.5 Units

The abbreviation of units for numeric results shall appear here. ISO standard abbreviations in accordance with ISO 2955⁴ should be employed when available (e.g., use mg rather than milligrams). Units can be reported in upper or lower case.

9.6 Reference Ranges

This value shall be reported in the following sample format: (lower limit to upper limit; i.e., 3.5 to 4.5). The range definition can be included by text description. See the following paragraph in this section. If a substance is toxic, then the upper limit of the range identifies the toxic limit. If the substance being measured is a drug, the lower limits identify the lower therapeutic bounds and the upper limits represent the upper therapeutic bounds above which toxic side effects are common.

A result may have multiple ranges, for example, an observation may have a physiological and a therapeutic range (e.g., serum magnesium is being used to treat eclampsia). When multiple ranges are sent, they shall be separated by repeat delimiters. Each range can also have a text description. The text description follows immediately after the range, and is separated from it by a component delimiter. Most results will only have one normal range transmitted.

9.7 Result Abnormal Flags

This field shall indicate the normalcy status of the result. The characters for representing significant changes either up or down or abnormal values shall be:

L	below low normal
H	above high normal
LL	below panic normal
HH	above panic high
<	below absolute low, that is off low scale on instrument
>	above absolute high, that is off high scale on an instrument
N	normal
A	abnormal
U	significant change up
D	significant change down
B	better, use when direction not relevant or not defined
W	worse, use when direction not relevant or not defined

When the instrument can discern the normal status of a textual report, such as microbiological culture, these should be reported as N when normal and A when abnormal.

9.8 Nature of Abnormality Testing

The kind of normal testing performed shall use the following representation: A denotes that an age-based population was tested; S, a sex-based population; and R, a race-based population. As many of the codes as apply shall be included. For example, if sex, age, and race normals were tested, an ASR would be transmitted. N implies that generic normal range was applied to all patient specimens.

9.9 Result Status

The following codes shall be used:

C	correction of previously transmitted results
P	preliminary results
F	final results
X	order cannot be done
I	in instrument, results pending
S	partial results
M	this result is an MIC level
R	this result was previously transmitted
N	this result record contains necessary information to run a new order

NOTE: For example, when ordering a sensitivity, the information system may download a result record containing the organism type, or species, identified in a previous test.

Q	this result is a response to an outstanding query
V	operator verified/approved result
W	Warning: Validity is questionable

9.10 Date of Change in Instrument Normative Values or Units

This field shall remain empty if there are no relevant normals or units. Otherwise, it shall be represented as in Section 5.6.2. A change in these data from those recorded in the receiving system's dictionary indicates a need for manual review of the results to detect whether they can be considered the same as preceding ones.

9.11 Operator Identification

The first component identifies the instrument operator who performed the test. The second component identifies the verifier for the test.

9.12 Date/Time Test Started

This field records the date and time the instrument started the test for which the results are now being reported. Date and time should be reported as specified in Section 5.6.2.

9.13 Date/Time Test Completed

This field records the date and time the instrument completed the test for which the results are now being reported. Date and time should be reported as specified in Section 5.6.2.

9.14 Instrument Identification

This field identifies the instrument or section of instrument that performed this particular measurement.

10 Comment Record

Comment records may be inserted anywhere except after the message terminator record. Each comment record shall apply to the first noncomment record preceding it. The comment record shall include the following fields:

10.1 Record Type ID

This record type shall be denoted by C.

10.2 Sequence Number

The sequence number is as defined in Section 5.6.7.

10.3 Comment Source

This field specifies the comment origination point as follows:

P	practice
L	information system
I	clinical instrument system

10.4 Comment Text

Where comment codes/mnemonics are used, the code should be sent first, followed, if desired, by the comment text and separated by a component delimiter as given in Section 5.6.6.

10.5 Comment Type

The following codes may be used to qualify comment record types:

G	generic/free result comment
T	result name comment
P	positive result comment
N	negative result comment
I	instrument flag(s) comment

11 Request Information Record

The request information record is used by either clinical instrument or information system to remotely request information from the reciprocal system.

NOTE: Only one request record may be outstanding at a time; the receiver of a request record must terminate the request, when finished, by means of the message terminator record, or the sender must cancel the request before sending a second logical request.

The request record shall include the following fields:

11.1 Record Type ID

This field is coded as Q.

11.2 Sequence Number

This field is as given in Section 5.6.7.

11.3 Starting Range ID Number

This field may contain three or more components to define a range of patients/specimens/manufacturers selection criteria. The first component is the information system patient ID number. The second component is the information system specimen ID number. Any further components are manufacturer-defined and for use in request subresult information (that is, an individual isolate/battery for a specimen

number). These components are position dependent. A list of sample IDs could be requested by the use of the repeat delimiter to separate IDs.

When ALL is entered, and the information system is sending the request record, it is taken to mean all specimen results ordered by the inquiring system. If the instrument is generating the request record, then it is taken to mean all demographics and tests being ordered should be sent to the instrument at this time. The request is then interpreted for that identified subset of specimens as further modified by the test specifications and date ranges as described below.

This specification does not address how long data are to be retained by an instrument, nor does it require that the instrument provide the search services implied by some of the field contents. The appropriate response for a request for results is simply the return of a subset of results that are currently in storage and can be practically retrieved.

11.4 Ending Range ID Number

This field is similar to that described in Section 11.3. If a single result or specimen demographic or test order is being requested, then this field may be left blank.

11.5 Universal Test ID

This field is as described in Section 5.6.1. This field may alternatively contain multiple codes separated by repeat delimiters, or the field may contain the text ALL, which signifies a request for all results on all tests or batteries for the patients/specimens/tests defined in Sections 11.3 and 11.4 and within the dates described in Sections 11.6 and 11.7.

11.6 Nature of Request Time Limits

Specify whether the date and time limits specified in Sections 11.7 and 11.8 refer to the specimen collect or ordered date (see Section 8.4.8) or test date (see Section 8.4.23): S indicates the specimen collect date; R indicates the result test date. If nothing is entered, the date criteria are assumed to be the result test date.

11.7 Beginning Request Results Date and Time

This field shall represent either a beginning (oldest) date and time for which results are being requested or a single date and time. The field may contain a single date and time or multiple dates and times separated by repeat delimiters. Each date and time shall be represented as specified in Section 5.6.2.

If no date and time is included, the instrument should assume that the information system wants results going as far into the past as possible and consistent with the criteria specified in other fields.

11.8 Ending Request Results Date and Time

This field, if not null, specifies the ending or latest (or most recent) date and time for which results are being requested. Date and time shall be represented as in Section 5.6.2.

11.9 Requesting Physician Name

This field identifies the individual physician requesting the results. The identity of the requesting physician is recorded as specified in Section 5.6.6.

11.10 Requesting Physician Telephone Number

This field is as specified in Section 5.6.3.

11.11 User Field Number 1

This is a user-defined field.

11.12 User Field Number 2

This is a user-defined field.

11.13 Request Information Status Codes

The following codes shall be used:

C	correction of previously transmitted results
P	preliminary results
F	final results
X	results cannot be done, request cancelled
I	request results pending
S	request partial/unfinalized results
M	result is an MIC level
R	this result was previously transmitted
A	abort/cancel last request criteria (allows a new request to follow)
N	requesting new or edited results only
O	requesting test orders and demographics only (no results)
D	requesting demographics only (e.g., patient record)

12 Message Terminator Record

This is the last record in the message. A header record may be transmitted after this record signifying the start of a second message.

12.1 Record Type ID

This record type is coded as L.

12.2 Sequence Number

This field is as described in Section 5.6.7. (For this record type, the value of this field should always be 1.)

12.3 Termination Code

This field provides an explanation of the end of the session.

Nil, N	normal termination
T	sender aborted
R	receiver requested abort
E	unknown system error
Q	error in last request for information

I no information available from last query
F last request for information processed

NOTE: I or Q will terminate a request and allow processing of a new request record.

13 Scientific Record

The scientific record exchanges the test data on clinical laboratory/instrument performance, quality assurance, or method development. It contains information in addition to the analyte measures found in the result record, although there are common elements in the two records.

13.1 Record Type ID

This field shall be identified by the character S.

13.2 Sequence Number

The sequence number shall be assigned as described in Section 5.6.7.

13.3 Analytical Method

This text field shall conform to Appendix I of Elevitch and Boroviczeny.⁵

13.4 Instrumentation

This text field shall be represented by an ID composed of the manufacturer and instrument codes connected by a dash (Latin-1(45)). These codes shall conform to Appendix I of Elevitch and Boroviczeny.⁵

13.5 Reagents

This text field shall include a list of constituent reagent codes, separated by subfield ID. These codes shall conform to the scheme of The American Chemical Society.

13.6 Units of Measure

The units of measure shall be represented as specified in Section 9.5.

13.7 Quality Control

Specifications for quality control are being developed.

13.8 Specimen Descriptor

This field shall use the convention described in Section 8.4.16.

13.9 Reserved Field

This field is reserved for future expansion.

13.10 Container

Specifications for containers are being developed.

13.11 Specimen ID

This text field shall represent a unique specimen identifier assigned by the originator and returned by the receiving instrument.

13.12 Analyte

Specifications for analyte are being developed.

13.13 Result

This numeric field shall represent the determined value of the analyte.

13.14 Result Units

This field shall be represented as described in Section 9.5.

13.15 Collection Date and Time

This field shall be represented in accordance with Section 5.6.2.

13.16 Result Date and Time

This field shall be represented in accordance with Section 5.6.2.

13.17 Analytical Preprocessing Steps

This text field shall contain the description of any preprocessing steps.

13.18 Patient Diagnosis

This field shall be represented as IDC-9-CM codes.

13.19 Patient Birthdate

This should be represented as specified in Section 7.8.

13.20 Patient Sex

This field shall be represented in accordance with Section 7.9.

13.21 Patient Race

This should be represented in accordance with Section 7.10.

14 Manufacturer Information Record

This record is provided solely for custom use by the instrument or computer system manufacturer. It has no inherent hierarchical level and may be inserted at any point except immediately following a message terminator record. It is recommended that this record type not be implemented unless all other possibilities have been exhausted. This record shall include the following:

- Record type ID—Coded as M.
- Sequence number—As defined in Section 5.6.7.

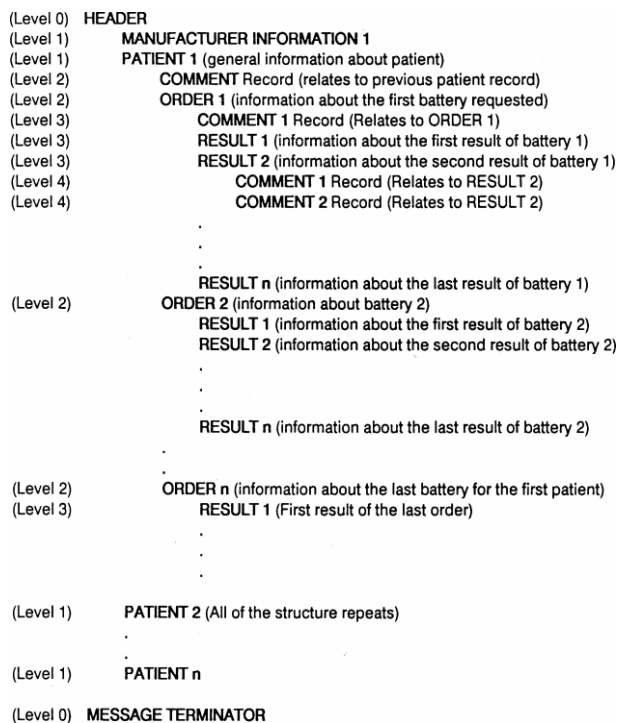


Figure 1. Logical Structure of a Message

Line #	Record Type	(Level) Increment	Action
A	Header	(Level 0)+0	
B	Patient 1	(Level 1)+1	
C	Order 1	(Level 2)+1	
D	Result 1	(Level 3)+1	
E	Order 2	(Level 2)-1	{ Save A-D }
F	Order 3	(Level 2)+0	
G	Patient 2	(Level 1)-1	{ Save E-F }
H	Order 1	(Level 2)+1	
I	Comment 1	(Level 3)+1	
J	Result 1	(Level 3)+0	
K	Comment 1	(Level 4)+1	
L	Result 2	(Level 3)-1	{ Save G-K }
M	Order 2	(Level 2)-1	{ Save L }
N	Patient 3	(Level 1)-1	{ Save M }
O	Order 1	(Level 2)+1	
P	Result 1	(Level 3)+1	
Q	Message Terminator	(Level 0)-3	{ Save N-P }

NOTE 1: Q is assumed as saved by virtue of the record type function.

NOTE 2: Given the following example, permanent storage of data, by the receiver, should occur at points: E, G, L, M, N, and Q.

Figure 2. Logical Information Storage Requirements

```

H|\^&<CR>
P|1<CR>
O|1|||^^^A1<CR>
R|1||0.356<CR>
P|2<CR>
O|1|||^^^A2<CR>
R|1||1.672<CR>
.
.
P|96<CR>
O|1|||^^^H12<CR>
R|1||0.402<CR>
L|1<CR>

```

NOTE 1: This sample is not recommended for implementation.

NOTE 2: Direction: instrument to information system.

Figure 3. Minimal Implementation (No Patient ID or Specimen ID)

```

H|\^&<CR>
P|1<CR>
O|1|927529|||^^^A1^^^A2<CR>
R|1|^^^A1|0.295|||||||19890327132247<CR>
R|2|^^^A2|0.312|||||||19890327132248<CR>
P|2<CR>
O|1|927533|||^^^A3^^^A4<CR>
R|1|^^^A3|1.121|||||||19890327132422<CR>
R|2|^^^A4|1.097|||||||19890317132423<CR>
L|1<CR>

```

Figure 4. No Patient ID; Specimen ID and Multiple Results Shown

```

H|\^&||PSWD|Harper Labs|2937 Southwestern Avenue^Buffalo^NY^73205||319 412-9722|||P|1394-97|19890314<CR>
O|1|^032989325|^032989327|ALL|||||||0<CR>
L|1|N<CR>

```

Figure 5. Request from Analyzer for Test Selections on Specimens 032989325-032989327

```

H|\^&||PSWD|Harper Labs|2937 Southwestern Avenue^Buffalo^NY^73205||319 412-9722|||P|1394-97|19890314<CR>
P|1|2734|123|306-87-4587|BLAKE^LINDSEY^ANN^MISS<CR>
O|1|032989325|||^^^BUN|R<CR>
O|2|032989325|||^^^ISE|R<CR>
O|3|032989325|||^^^HDL^^^GLU|R<CR>
P|2|2462|158|287-17-2791|POHL^ALLEN^M.^<CR>
O|1|032989326|||^^^LIVER^^^GLU|S<CR>
P|3|1583|250|151-37-6926|SIMPSON^ALBERT^^MR<CR>
O|1|032989327|||^^^CHEM12^^^LIVER|R<CR>
L|1|F<CR>

```

Figure 6. Response from Information System for Previous Request

```

H|\&||PSWD|Harper Labs|2937 Southwestern Avenue^Buffalo^NY^73205||319 412-9722|||P|1394-97|19890314<CR>
P|1|2734|123|306-87-4587|BLAKE^LINDSEY^ANN^MISS<CR>
C|1|L|Notify IDC if tests positive|G<CR>
O|1|032989325|^^^BUN|R<CR>
R|1|^^^BUN|8.71<CR>
C|1|I|TGP^Test Growth Positive|P<CR>
C|2|I|colony count >10,000|P<CR>
O|2|032989325|^^^ISE|R<CR>
R|1|^^^ISE^NA|139|mEq/L<CR>
R|2|^^^ISE^K|4.2|mEq/L<CR>
R|3|^^^ISE^CL|111|mEq/L<CR>

O|3|032989325|^^^HDL|R<CR>
R|1|^^^HDL|70.29<CR>
O|4|032989325|^^^GLU|R<CR>
R|1|^^^GLU|92.98<CR>
C|1|I|Reading is Suspect|I<CR>
P|2|2462|158|287-17-2791|POHL^ALLEN^M.<CR>
O|1|032989326|^^^LIVER|S<CR>
R|1|^^^LIVER^AST|29<CR>
R|2|^^^LIVER^ALT|50<CR>
R|3|^^^LIVER^TBILI|7.9<CR>
R|4|^^^LIVER^GGT|29<CR>

O|2|032989326|^^^GLU|S<CR>
R|1|^^^GLU|91.5<CR>
P|3|1583|250|151-37-6926|SIMPSON^ALBERT^MR<CR>
O|1|032989327|LIVER|R<CR>
R|1|^^^AST|28<CR> [Test ID field Implicitly Relates to LIVER order 1]
R|2|^^^ALT|49<CR>
R|3|^^^TBILI|7.3<CR>
R|4|^^^GGT|27<CR>

O|2|032989327|CHEM12|R<CR>
R|1|^^^CHEM12^ALB-G|28<CR> [Test ID field Explicitly Relates to CHEM12 order 1]
R|2|^^^CHEM12^BUN|49<CR>
R|3|^^^CHEM12^CA|7.3<CR>
R|4|^^^CHEM12^CHOL|27<CR>
R|5|^^^CHEM12^CREAT|4.2<CR>
R|6|^^^CHEM12^PHOS|12<CR>
R|7|^^^CHEM12^GLUHK|9.7<CR>
R|8|^^^CHEM12^NA|138.7<CR>
R|9|^^^CHEM12^K|111.3<CR>
R|10|^^^CHEM12^CL|6.7<CR>
R|11|^^^CHEM12^UA|7.3<CR>
R|12|^^^CHEM12^TP|9.2<CR>

L|1<CR>

```

Figure 7. Results from Given Ordered Test Selections Shown in Various Formats

```

H|\&||PSWD|Harper Labs|2937 Southwestern Avenue^Buffalo^NY^73205||319 412-9722|||P|1394-97|19890314<CR>
Q|1|032989326|ALL<CR>
L|1<CR>

```

Figure 8. Request from Information System to Instrument for Previously Run Results

```

H|^&||PSWD|Harper Labs|2937 Southwestern Avenue^Buffalo^NY^73205||319 412-9722|||P|1394-97|19890314<CR>
P|1|2462|158|287-17-2791|POHL^ALLEN^M.<CR>
O|1|032989326||^LIVER|S<CR>
R|1|^AST|29<CR>
R|2|^ALT|50<CR>
R|3|^TBILI|7.9<CR>
R|4|^GGT|29<CR>
O|2|032989326||^GLU|S<CR>
R|1|^GLU|91.5<CR>
L|1|F<CR>

```

Figure 9. Reply to Result Request

```

H|^&||Password1|Micro1||||LSI1||P|1394-94|19890501074500<CR>
P|1||52483291||Smith^John|Samuels|19600401|M|W|4526 C Street^Fresno^CA^92304|||(402)782-
3424X242|542^Dr. Brown||72^in.|175^lb.||Penicillin||||19890428|IP|Ward1|||C|M|WSP||ER|PC^Prompt Care<CR>
O|1|5762^01||^BC^Blood Culture^POSCOMBO|R|198905011530|198905020700||456^Farnsworth|N||-
198905021130|BL^Blood|123^Dr. Wirth|||||Instrument#1|||ER|N<CR>
R|1|^Org#|51^Strep Species||N<CR>
R|2|^Bio|BH+^Beta Hemolytic||N<CR>
L|1<CR>

```

Figure 10. Microbiology Order and Result-Download of Demographics and Order

```

H|^&||Password1|Micro1||||LSI1||P|1394-94|19890501074500<CR>
P|1||52483291<CR>
O|1|5762^01||^BC^POSCOMBO|||||||BL|||||||F<CR>
R|1|^ORG#|103^Group D Enteroc<CR>
R|2|^AM^MIC|>16<CR>
R|3|^AM^INTERP1|++<CR>
R|4|^AM^DOSAGE1|PO 250-500 mg Q6h<CR>
R|5|^AM^DOSAGE1^COSTCODE|$25<CR>
R|6|^AM^INTERP2|+++<CR>
R|7|^AM^DOSAGE2|IV 1.0-2.0 gm Q4h<CR>
R|8|^P^MIC|<0.25<CR>
R|9|^P^INTERP1|++<CR>
R|10|^P^DOSAGE1|PO 250-500 mg Q6h<CR>
R|11|^P^DOSAGE1^COSTCODE|$25<CR>
R|12|^P^INTERP2|+++<CR>
R|13|^P^DOSAGE2|IM 0.9-1.2 MIL U Q6-12h<CR>
.
.
.
R|90|^BIOTYPE|102-34021<CR>
L|1<CR>

```

Figure 11. Microbiology Order and Result-Upload of Finalized Results

References

- ¹ ISO. *Information Processing—8-bit single-byte coded graphic character sets—Part 1: Latin Alphabet No. 1*. ISO 8859-1. Geneva: International Organization of Standardization; 1987.
- ² ANSI X3.30 *Representation for Calendar Date and Ordinal Date for Information Interchange*. New York, NY: American National Standards Institute.
- ³ ANSI X3.43 *Representations of Local Time of Day for Information Interchange*. New York, NY: American National Standards Institute.
- ⁴ ISO. *Information Processing—Representation of SI and Other Units in Systems with Limited Character Sets*. ISO 2955. Geneva: International Organization for Standardization; 1993.
- ⁵ Elevitch FR, Boroviczeny KG. *Transfer Data: A Proposed International Standard for Interlaboratory Information Exchange*. Skokie, IL: College of American Pathologists.

Appendix. Comparison of LIS02 and LIS5

A1. Table A1 shows the major differences in requirements between NCCLS documents LIS02—*Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems* and LIS5—*Standard Specification for Transferring Clinical Observations Between Independent Computer Systems*. Other modifications and additions have been made. Not all of the fields required in NCCLS document LIS5 are required in this specification. It is the responsibility of the user of this standard to compare the requirements of these two specifications and the changes that have been incorporated since the last issue.

Table A1. Requirement Comparison Between LIS02 and LIS5

Requirement	LIS02	LIS5
Addenda record	No	Yes (see Sections 5.2 and 6.5.9)
Maximum record length	No	Yes (see Section 6.2 and Table 1)
Maximum field length	No	Yes (see Section 6.2 and Table 1)
Error check record	No	Yes (see Section 5.3)
Extended character set	Yes (see Section 5.1)	No
Latin-1(10) linefeed	No	Yes (see Section 6.1.1)
Record Type M	Yes (see Section 5.4)	No
Repeat delimiter	Yes (see Section 5.4.3)	*
Component delimiter	Yes (see Section 5.4.4)	†
Escape delimiter	Yes (see Section 5.4.5)	No
Linked results	No	Yes (see Section 9.4.27)
Multiple specimen orders	No	Yes (see Section 9.4.29)
Universal Test ID	4-part‡ (see Section 5.6.1)	3-part (see Section 6.6.1)
Time fields	HHMMSS with optional time zone (see Section 5.6.2)	HHMM (see Section 6.6.2.1)
Result record representation§	Single analytic result value (see Section 9.4)	Multiple parts permitted (see Section 10.1.4.2)
Manufacturer's or local code	Yes—must be defined and documented by manufacturer (see Section 5.6.1.4)	No
Source of specimen	Specimen descriptor (see Section 8.4.16)	Specimen type and location (see Section 9.4.16)

* In NCCLS document LIS5, this delimiter is referred to as a "subfield delimiter."

† In NCCLS document LIS5, this delimiter is referred to as a "sub-subfield delimiter."

‡ Manufacturer-defined test code added.

§ This record was modified and renamed in LIS02.

NCCLS consensus procedures include an appeals process that is described in detail in Section 8 of the Administrative Procedures. For further information contact the Executive Offices or visit our website at www.nccls.org.

Summary of Delegate Comments and Area Committee Responses

LIS02-A2: *Specification for Transferring Information Between Clinical Laboratory Instruments and Information Systems; Approved Standard—Second Edition*

General

1. I have some comments on the entire concept. It seems to be very optimistic and ambitious to expect this much data being transferred between instrument and LIS (after all the result reports are coming out of the LIS and not the instrument). To have the ability to transfer all information, even down to height of patient in centimeters, will be extremely costly and time consuming for LIS vendors, LIS administrators, and vendors of the instruments. Unless the instrument is acting as the LIS repository for data, and the result reports are being generated by the instrument, I believe this to be overkill. The experience I have had with interfacing instruments has included several identifiers such as patient name, specimen number, and date of birth. This should be adequate to match the results coming off of the instrument to the record containing the rest of the demographics, etc. in the LIS.

The “HL7 format” that NCCLS seems to have modeled this document over is a logical approach for communicating between instrument and LIS. The fields in the message need to be limited, as mentioned previously, however. I have not had time to compare HL7 to the format that was described in the document, but they need to make sure that if they are going to use this type of format, that it is the same as the HL7 format for reporting and interoperability that PHIN compliance is pushing. It will be much easier for vendors of commercial of the shelf (COTS) systems to incorporate this functionality into their systems if they are already using the identical format. A different format will require more developmental time and costs.

One other note on the amount of data fields they are suggesting: Most instrumentation does not require specifics on diagnosis, height, weight, race, etc. to be inputted for a result to be calculated. Actually, I would guess that this is a rarity if even applicable. I would not want to have NCCLS suggesting that all instrument interfaces be set up to this extent, when it is not needed. To suggest a format if it is needed is okay, but only if it is needed.

To sum it up: I feel the HL7 format is the logical choice so that there is a “standard” the vendors can utilize. The extent of this needs to be limited by application, otherwise the concept will fail due to the high demand on resources and the costs involved.

- **Response to Comment 1, paragraphs 1 and 3:**

The majority of the fields are optional (i.e., it is not necessary for the sender to send this information and the receiver may ignore this information if it is not relevant for the processing). The existing implementations of the ASTM 1394 standard (predecessor of LIS02) utilize these mechanisms. However, some instruments may need the data (e.g., the mentioned patient weight, height) for determination of proper ranges in the validation of results or other purposes. The area committee believes that it is correct to leave the potential capabilities and that the optional use will not lead to “overkill” as mentioned in the comment.

Response to Comment 1, paragraph 2:

Historically, the message format of HL7 and ASTM 1394 standards has the same origin in the 1980s. However, in the meantime, these standards were developed independently and the scope for both standards is currently very different. For the LIS02 document, as a successor of the ASTM 1394 standard, it is more important and relevant to keep the backwards compatibility with this standard instead of trying to follow HL7 development.

Response to Comment 1, paragraph 4: NCCLS document AUTO3—*Laboratory Automation: Communications with Automated Clinical Laboratory Systems, Instruments, Devices, and Information Systems* uses HL7 as a base and keeps the compatibility with the HL7 format.

2. Specifications seem more appropriate to a data manager system that functions as a mini-LIS. The instruments interfaced to our LIS do not transmit most of the fields specified.

- **The majority of the fields are optional (i.e., it is not necessary for the sender to send this information and the receiver may ignore this information if it is not relevant for the processing). The existing implementations of the ASTM 1394 standard (predecessor of LIS02) utilize these mechanisms. However, some instruments may need the data (e.g., the mentioned patient weight, height) for determination of proper ranges in the validation of results or other purposes. We feel it is correct to leave the potential capabilities for optional use in this document.**

Section 5.4, Delimiters

3. “Repeat delimiter = backslash () Latin-1 (96)” should be: Repeat delimiter = backslash (\) Latin-1 (96).

- **The typographical error has been corrected.**

Section 13.12, Analyte

4. The document states “Specifications for analyte are being developed.” Will this be completed before the document has its final acceptance?
- **The development of vocabulary specifications for analyte is beyond the scope of this standard.**

The Quality System Approach

NCCLS subscribes to a quality system approach in the development of standards and guidelines, which facilitates project management; defines a document structure via a template; and provides a process to identify needed documents. The approach is based on the model presented in the most current edition of NCCLS document HS1—*A Quality System Model for Health Care*. The quality system approach applies a core set of “quality system essentials” (QSEs), basic to any organization, to all operations in any healthcare service’s path of workflow (i.e., operational aspects that define how a particular product or service is provided). The QSEs provide the framework for delivery of any type of product or service, serving as a manager’s guide. The quality system essentials (QSEs) are:

Documents & Records	Equipment	Information Management	Process Improvement
Organization	Purchasing & Inventory	Occurrence Management	Service & Satisfaction
Personnel	Process Control	Assessment	Facilities & Safety

LIS02-A2 addresses the quality system essentials (QSEs) indicated by an “X.” For a description of the other NCCLS documents listed in the grid, please refer to the Related NCCLS Publications section on the following page.

Documents & Records	Organization	Personnel	Equipment	Purchasing & Inventory	Process Control	Information Management	Occurrence Management	Assessment	Process Improvement	Service & Satisfaction	Facilities & Safety
			LIS3	LIS3	X LIS4	X AUTO3 LIS1 LIS5 LIS8	LIS6				

Adapted from NCCLS document HS1—*A Quality System Model for Health Care*.

Path of Workflow

A path of workflow is the description of the necessary steps to deliver the particular product or service that the organization or entity provides. For example, GP26-A2 defines a clinical laboratory path of workflow which consists of three sequential processes: preanalytic, analytic, and postanalytic. All clinical laboratories follow these processes to deliver the laboratory’s services, namely quality laboratory information.

LIS02-A2 addresses the clinical laboratory path of workflow steps indicated by an “X.” For a description of the other NCCLS documents listed in the grid, please refer to the Related NCCLS Publications section on the following page.

Preanalytic					Analytic		Postanalytic	
Patient Assessment	Test Request	Specimen Collection	Specimen Transport	Specimen Receipt	Testing Review	Laboratory Interpretation	Results Report	Post-test Specimen Management
							X	

Adapted from NCCLS document HS1—*A Quality System Model for Health Care*.

Related NCCLS Publications

- AUTO3-A** **Laboratory Automation: Communications with Automated Clinical Laboratory Systems, Instruments, Devices, and Information Systems; Approved Standard (2000).** This document provides standards to facilitate accurate and timely electronic exchange of data and information between the automated laboratory elements.
- LIS1-A** **Standard Specification for Low-Level Protocol to Transfer Messages Between Clinical Laboratory Instruments and Computer Systems (2003).** This specification describes the electronic transmission of digital information between the clinical laboratory instruments (those that measure one or more parameters from one or multiple samples) and computer systems (those that are configured to accept instrument results for further processing, storage, reporting, or manipulation).
- LIS3-A** **Standard Guide for Selection of a Clinical Laboratory Information Management System (2003).** This guide covers the selection, purchase, use, enhancement, and updating of computer technology supplied by a vendor as a complete system in the clinical laboratory. The purpose of the guide is to assist hospitals, clinics, and independent laboratories through the entire automation project in order to minimize the risks and maximize the benefits. It also includes checklists of items and design aids to be considered at each stage of planning to assist in carrying out the project.
- LIS4-A** **Standard Guide for Documentation of Clinical Laboratory Computer Systems (2003).** This guide covers documentation (defined as the information needed to install, use, maintain, or modify the system) for a computer system operating in a clinical laboratory.
- LIS5-A** **Standard Specification for Transferring Clinical Observations Between Independent Computer Systems (2003).** This specification details how clinical observations can be transferred between independent computer systems.
- LIS6-A** **Standard Practice for Reporting Reliability of Clinical Laboratory Information Systems (2003).** This practice describes a system for collecting data, maintaining records, and reporting on the reliability of operating clinical laboratory computer systems. The reliability measure will be achieved by documenting the number, severity, cause, impact, and duration of the failures that a system experiences. This practice can be implemented with paper forms or computer records.
- LIS8-A** **Standard Guide for Functional Requirements of Clinical Laboratory Information Management Systems (2003).** This guide covers the capabilities needed for a Clinical Laboratory Information Management System (CLIMS). It was written so that both vendors/developers of CLIMS and laboratory managers would have a common understanding of the requirements and logical structure of a laboratory data system. This guide will also provide more uniformity in the way that requirements are expressed from one laboratory to another.

NOTES

Active Membership (as of 1 October 2004)

Sustaining Members

Abbott Laboratories
American Association for Clinical Chemistry
Bayer Corporation
BD
Beckman Coulter, Inc.
bioMérieux, Inc.
CLMA
College of American Pathologists
GlaxoSmithKline
Ortho-Clinical Diagnostics, Inc.
Pfizer Inc
Roche Diagnostics, Inc.

Professional Members

American Academy of Family Physicians
American Association for Clinical Chemistry
American Association for Respiratory Care
American Chemical Society
American Medical Technologists
American Society for Clinical Laboratory Science
American Society for Microbiology
American Society of Hematology
American Type Culture Collection, Inc.
Asociación Mexicana de Bioquímica Clínica A.C.
Assn. of Public Health Laboratories
Assoc. Micro. Clinici Italiani-A.M.C.L.I.
British Society for Antimicrobial Chemotherapy
Canadian Society for Medical Laboratory Science - Société Canadienne de Science de Laboratoire Médical
Canadian Standards Association
Clinical Laboratory Management Association
COLA
College of American Pathologists
College of Medical Laboratory Technologists of Ontario
College of Physicians and Surgeons of Saskatchewan
ESCMID
International Council for Standardization in Haematology
International Federation of Biomedical Laboratory Science
International Federation of Clinical Chemistry
Italian Society of Clinical Biochemistry and Clinical Molecular Biology
Japan Society of Clinical Chemistry
Japanese Committee for Clinical Laboratory Standards
Joint Commission on Accreditation of Healthcare Organizations
National Academy of Clinical Biochemistry
National Association of Testing Authorities - Australia
National Society for Histotechnology, Inc.
New Zealand Association of Phlebotomy
Ontario Medical Association Quality Management Program-Laboratory Service
RCPA Quality Assurance Programs
PTY Limited
Sociedad Española de Bioquímica Clínica y Patología Molecular
Sociedade Brasileira de Análises Clínicas
Taiwanese Committee for Clinical Laboratory Standards (TCCLS)
Turkish Society of Microbiology

Government Members

Armed Forces Institute of Pathology
Association of Public Health Laboratories
BC Centre for Disease Control
Caribbean Epidemiology Centre
Centers for Disease Control and Prevention
Centers for Medicare & Medicaid Services
Centers for Medicare & Medicaid Services/CLIA Program

Chinese Committee for Clinical Laboratory Standards
Commonwealth of Pennsylvania Bureau of Laboratories
Department of Veterans Affairs
Deutsches Institut für Normung (DIN)
FDA Center for Devices and Radiological Health
FDA Center for Veterinary Medicine
FDA Division of Anti-Infective Drug Products
Iowa State Hygienic Laboratory
Massachusetts Department of Public Health Laboratories
National Center of Infectious and Parasitic Diseases (Bulgaria)
National Health Laboratory Service (South Africa)
National Institute of Standards and Technology
National Pathology Accreditation Advisory Council (Australia)
New York State Department of Health
Ontario Ministry of Health
Pennsylvania Dept. of Health
Saskatchewan Health-Provincial Laboratory
Scientific Institute of Public Health; Belgium Ministry of Social Affairs, Public Health and the Environment
Swedish Institute for Infectious Disease Control

Industry Members

AB Biodisk
Abbott Laboratories
Abbott Diabetes Care
Acrometrix Corporation
Advancis Pharmaceutical Corporation
Affymetrix, Inc.
Ammirati Regulatory Consulting
Anna Longwell, PC
A/S ROSCO
AstraZeneca Pharmaceuticals
Aventis
Axis-Shield POC AS
Bayer Corporation - Elkhart, IN
Bayer Corporation - Tarrytown, NY
Bayer Corporation - West Haven, CT
BD
BD Consumer Products
BD Diagnostic Systems
BD Thailand Ltd.
BD VACUTAINER Systems
Beckman Coulter, Inc.
Beckman Coulter K.K. (Japan)
Bio-Development SRL
Bio-Inova Life Sciences International
Biomedica Laboratories SDN BHD
bioMérieux, Inc. (MO)
Biometrology Consultants
Bio-Rad Laboratories, Inc.
Bio-Rad Laboratories, Inc. – France
Bio-Rad Laboratories, Inc. – Plano, TX
Blaine Healthcare Associates, Inc.
Bristol-Myers Squibb Company
Canadian External Quality Assessment Laboratory
Cepheid
Chen & Chen, LLC
Chiron Corporation
ChromaVision Medical Systems, Inc.
Clinical Micro Sensors
The Clinical Microbiology Institute
Cognigen
CONOSCO
Copan Diagnostics Inc.
Cosmetic Ingredient Review
Cubist Pharmaceuticals
Dade Behring Inc. - Cupertino, CA
Dade Behring Inc. - Deerfield, IL
Dade Behring Inc. - Glasgow, DE
Dade Behring Inc. - Marburg, Germany
Dade Behring Inc. - Sacramento, CA
David G. Rhoads Associates, Inc.
Diagnostic Products Corporation
Digene Corporation
Eiken Chemical Company, Ltd.
Elanco Animal Health
Electa Lab s.r.l.
Enterprise Analysis Corporation
F. Hoffman-La Roche AG
Fort Dodge Animal Health
Gen-Probe
GenVault
GlaxoSmithKline
Greiner Bio-One Inc.
Immunic Corporation
ImmunoSite, Inc.
Instrumentation Laboratory
International Technidyne Corporation
I-STAT Corporation
Johnson and Johnson Pharmaceutical Research and Development, L.L.C.
K.C.J. Enterprises
LifeScan, Inc. (a Johnson & Johnson Company)
Machao Diagnostics
Medical Device Consultants, Inc.
Merck & Company, Inc.
Minigrip/Zip-Pak
Nanosphere, Inc.
National Pathology Accreditation Advisory Council (Australia)
Nippon Becton Dickinson Co., Ltd.
Nissui Pharmaceutical Co., Ltd.
Norfolk Associates, Inc.
Novartis Pharmaceuticals Corporation
Olympus America, Inc.
Optimer Pharmaceuticals, Inc.
Ortho-Clinical Diagnostics, Inc. (Rochester, NY)
Ortho-McNeil Pharmaceutical (Raritan, NJ)
Oxoid Inc.
Paratek Pharmaceuticals
Pfizer Animal Health
Pfizer Inc
Pfizer Italia Srl
Powers Consulting Services
Premier Inc.
Procter & Gamble Pharmaceuticals, Inc.
QSE Consulting
Radiometer America, Inc.
Radiometer Medical A/S
Replidyne
Roche Diagnostics GmbH
Roche Diagnostics, Inc.
Roche Laboratories (Div. Hoffmann-La Roche Inc.)
Sarstedt, Inc.
Schering Corporation
Schleicher & Schuell, Inc.
Second Opinion
SFBC Anapharm
Streck Laboratories, Inc.
SYN X Pharma Inc.
Sysmex Corporation (Japan)
Sysmex Corporation (Long Grove, IL)
TheraDoc
Theravance Inc.
Thrombodyne, Inc.
THYMED GmbH
Transasia Engineers
Trek Diagnostic Systems, Inc.
Vetequinol S.A.
Vicuron Pharmaceuticals Inc.
Vysis, Inc.
Wyeth Research
XDX, Inc.
YD Consultant
YD Diagnostics (Seoul, Korea)

Trade Associations

AdvaMed
Japan Association of Clinical Reagents Industries (Tokyo, Japan)

Associate Active Members

Academisch Ziekenhuis - VUB (Belgium)
Alfred I. duPont Hospital for Children (DE)
All Children's Hospital (FL)
Allina Health System (MN)
American University of Beirut Medical Center (NY)
Anne Arundel Medical Center (MD)
Antwerp University Hospital (Belgium)
Arkansas Department of Health
ARUP at University Hospital (UT)
Associated Regional & University Pathologists (UT)
Atlantic Health System (NJ)
Aurora Consolidated Laboratories (WI)
AZ Sint-Jan (Belgium)

Azienda Ospedale Di Lecco (Italy)
Barnes-Jewish Hospital (MO)
Baxter Regional Medical Center (AR)
Baystate Medical Center (MA)
Bhagwas Duzen Laboratories (Turkey)
BC Biomedical Laboratories (Surrey, BC, Canada)
Bermuda Hospitals Board
Brazosport Memorial Hospital (TX)
Broward General Medical Center (FL)
Cadham Provincial Laboratory (Winnipeg, MB, Canada)
Calgary Laboratory Services (Calgary, AB, Canada)
California Pacific Medical Center
Canterbury Health Laboratories (New Zealand)
Cape Breton Healthcare Complex (Nova Scotia, Canada)
Carlilion Consolidated Laboratory (VA)
Carolinas Medical Center (NC)
Cathay General Hospital (Taiwan)
Central Texas Veterans Health Care System
Centro Diagnostico Italiano (Milano, Italy)
Champlain Valley Physicians Hospital (NY)
Chang Gung Memorial Hospital (Taiwan)
Changi General Hospital (Singapore)
Children's Hospital (NE)
Children's Hospital & Clinics (MN)
Children's Hospital Medical Center (Akron, OH)
Children's Medical Center of Dallas (TX)
Chinese Association of Advanced Blood Bankers (Beijing)
CHR St. Joseph Warquignies (Belgium)
Clarian Health - Methodist Hospital (IN)
CLSI Laboratories (PA)
Community Hospital of Lancaster (PA)
Community Hospital of the Monterey Peninsula (CA)
CompuNet Clinical Laboratories (OH)
Cook Children's Medical Center (TX)
Cook County Hospital (IL)
Covance Central Laboratory Services (IN)
Creighton University Medical Center (NE)
Danish Veterinary Laboratory (Denmark)
Detroit Health Department (MI)
DFS/CLIA Certification (NC)
Diagnósticos da América S/A (Brazil)
Dr. Everett Chalmers Hospital (New Brunswick, Canada)
Duke University Medical Center (NC)
Dwight David Eisenhower Army Medical Center (GA)
Eastern Health Pathology (Australia)
Emory University Hospital (GA)
Enzo Clinical Labs (NY)
Evangelical Community Hospital (PA)
Fairview-University Medical Center (MN)
Florida Hospital East Orlando
Focus Technologies (CA)
Focus Technologies (VA)
Foothills Hospital (Calgary, AB, Canada)
Franciscan Shared Laboratory (WI)
Fresno Community Hospital and Medical Center
Gamma Dynacare Medical Laboratories (Ontario, Canada)
Geisinger Medical Center (PA)
Guthrie Clinic Laboratories (PA)
Hagerstown Medical Laboratory (MD)
Harris Methodist Fort Worth (TX)
Hartford Hospital (CT)
Headwaters Health Authority (Alberta, Canada)
Health Network Lab (PA)
Health Partners Laboratories (VA)

Highlands Regional Medical Center (FL)	Maimonides Medical Center (NY)	Providence Health Care (Vancouver, BC, Canada)	Swedish Medical Center - Providence Campus (WA)
Hoag Memorial Hospital Presbyterian (CA)	Marion County Health Department (IN)	Provincial Laboratory for Public Health (Edmonton, AB, Canada)	Temple University Hospital (PA)
Holy Cross Hospital (MD)	Martin Luther King/Drew Medical Center (CA)	Queen Elizabeth Hospital (Prince Edward Island, Canada)	Tenet Odessa Regional Hospital (TX)
Hôpital du Sacré-Coeur de Montreal (Montreal, Quebec, Canada)	Massachusetts General Hospital (Microbiology Laboratory)	Queensland Health Pathology Services (Australia)	The Toledo Hospital (OH)
Hôpital Maisonneuve - Rosemont (Montreal, Canada)	MDS Metro Laboratory Services (Burnaby, BC, Canada)	Quest Diagnostics Incorporated (CA)	Touro Infirmary (LA)
Hôpital Saint-Luc (Montreal, Quebec, Canada)	Medical College of Virginia Hospital	Quintiles Laboratories, Ltd. (GA)	Tripler Army Medical Center (HI)
Hospital Consolidated Laboratories (MI)	Medical University of South Carolina	Regions Hospital	Truman Medical Center (MO)
Hospital for Sick Children (Toronto, ON, Canada)	Memorial Medical Center (Napoleon Avenue, New Orleans, LA)	Rex Healthcare (NC)	UCLA Medical Center (CA)
Hospital de Sousa Martins (Portugal)	Memorial Medical Center (Jefferson Davis Pkwy., New Orleans, LA)	Rhode Island Department of Health Laboratories	UCSF Medical Center (CA)
Hotel Dieu Grace Hospital (Windsor, ON, Canada)	Methodist Hospital (Houston, TX)	Robert Wood Johnson University Hospital (NJ)	UNC Hospitals (NC)
Huddinge University Hospital (Sweden)	Methodist Hospital (San Antonio, TX)	Royal Columbian Hospital (New Westminster, BC, Canada)	Unidad de Patología Clínica (Mexico)
Hunter Area Health Service (Australia)	Michigan Department of Community Health	Sahlgrenska Universitetssjukhuset (Sweden)	Union Clinical Laboratory (Taiwan)
Indiana University	Mid America Clinical Laboratories, LLC (IN)	St. Alexius Medical Center (ND)	Universita Campus Bio-Medico (Italy)
Innova Fairfax Hospital (VA)	Middlesex Hospital (CT)	St. Anthony Hospital (CO)	University College Hospital (Galway, Ireland)
Institute of Medical and Veterinary Science (Australia)	Monmouth Medical Center (NJ)	St. Anthony's Hospital (FL)	University of Alabama-Birmingham Hospital
International Health Management Associates, Inc. (IL)	Montreal Children's Hospital (Canada)	St. Barnabas Medical Center (NJ)	University of Chicago Hospitals (IL)
Jackson Memorial Hospital (FL)	Montreal General Hospital (Canada)	St. Christopher's Hospital for Children (PA)	University of Colorado Hospital
Jacobi Medical Center (NY)	National Serology Reference Laboratory (Australia)	St-Eustache Hospital (Quebec, Canada)	University of Debrecen Medical Health and Science Center (Hungary)
John C. Lincoln Hospital (AZ)	NB Department of Health & Wellness (New Brunswick, Canada)	St. John Hospital and Medical Center (MI)	University of Illinois Medical Center
Johns Hopkins Medical Institutions (MD)	The Nebraska Medical Center	St. John's Hospital & Health Center (CA)	University of Maryland Medical System
Kadlec Medical Center (WA)	New Britain General Hospital (CT)	St. Joseph Mercy Hospital (MI)	University of the Ryukyus (Japan)
Kaiser Permanente (MD)	New England Fertility Institute (CT)	St. Joseph's Hospital – Marshfield Clinic (WI)	University of Wisconsin Hospital
Kangnam St. Mary's Hospital (Korea)	New England Medical Center (MA)	St. Jude Children's Research Hospital (TN)	University of Washington
Kantonsspital (Switzerland)	New York City Department of Health & Mental Hygiene	St. Mary of the Plains Hospital (TX)	USA MEDDAC-AK
Kenora-Rainy River Regional Laboratory Program (Ontario, Canada)	NorDx (ME)	St. Michael's Hospital (Toronto, ON, Canada)	US LABS, Inc. (CA)
Kimball Medical Center (NJ)	North Carolina State Laboratory of Public Health	Ste. Justine Hospital (Montreal, PQ, Canada)	UZ-KUL Medical Center (Belgium)
King Abdulaziz Medical City – Jeddah (Saudi Arabia)	North Central Medical Center (TX)	Salem Clinic (OR)	VA (Hampton) Medical Center (VA)
King Faisal Specialist Hospital (Saudi Arabia)	North Shore - Long Island Jewish Health System Laboratories (NY)	San Francisco General Hospital (CA)	VA (Tuskegee) Medical Center (AL)
LabCorp (NC)	North Shore University Hospital (NY)	Santa Clara Valley Medical Center (CA)	Valley Children's Hospital (CA)
Laboratoire de Santé Publique du Quebec (Canada)	Northwestern Memorial Hospital (IL)	Seoul Nat'l University Hospital (Korea)	Virginia Beach General Hospital (VA)
Laboratorio Dr. Echevarne (Spain)	Ochsner Clinic Foundation (LA)	Shands at the University of Florida	Virginia Department of Health
Laboratório Fleury S/C Ltda. (Brazil)	O.L. Vrouwziekenhuis (Belgium)	South Bend Medical Foundation (IN)	Virginia Regional Medical Center (MN)
Laboratorio Manlab (Argentina)	Ordre professionnel des technologistes médicaux du Québec	South Western Area Pathology Service (Australia)	ViroMed Laboratories (MN)
Laboratory Corporation of America (NJ)	Orlando Regional Healthcare System (FL)	Southern Maine Medical Center	Washington Adventist Hospital (MD)
LAC and USC Healthcare Network (CA)	Ospedali Riuniti (Italy)	Spartanburg Regional Medical Center (SC)	Washoe Medical Center Laboratory (NV)
Lakeland Regional Medical Center (FL)	The Ottawa Hospital (Ottawa, ON, Canada)	Specialty Laboratories, Inc. (CA)	Waterford Regional Hospital (Ireland)
Landstuhl Regional Medical Center (APO AE)	OU Medical Center (OK)	State of Connecticut Dept. of Public Health	Wellstar Health Systems (GA)
LeBonheur Children's Medical Center (TN)	Our Lady of the Resurrection Medical Center (IL)	State of Washington Department of Health	West Jefferson Medical Center (LA)
Lewis-Gale Medical Center (VA)	Pathology and Cytology Laboratories, Inc. (KY)	Stony Brook University Hospital (NY)	Wilford Hall Medical Center (TX)
L'Hotel-Dieu de Quebec (Canada)	Pathology Associates Medical Laboratories (WA)	Stormont-Vail Regional Medical Center (KS)	William Beaumont Army Medical Center (TX)
Libero Instituto Univ. Campus BioMedico (Italy)	Peking University Shenzhen Hospital (China)	Sun Health-Boswell Hospital (AZ)	William Beaumont Hospital (MI)
Loma Linda Mercantile (CA)	The Permanente Medical Group (CA)	Sunnybrook Health Science Center (ON, Canada)	William Osler Health Centre (Brampton, ON, Canada)
Long Beach Memorial Medical Center (CA)	Piedmont Hospital (GA)	Sunrise Hospital and Medical Center (NV)	Winn Army Community Hospital (GA)
Louisiana State University Medical Center	Pocono Medical Center (PA)		Winnipeg Regional Health Authority (Winnipeg, Canada)
Lourdes Hospital (KY)			Wishard Memorial Hospital (IN)

OFFICERS

Thomas L. Hearn, Ph.D.,
President
Centers for Disease Control and Prevention

Robert L. Habig, Ph.D.,
President Elect
Abbott Laboratories

Wayne Brinster,
Secretary
BD

Gerald A. Hoeltge, M.D.,
Treasurer
The Cleveland Clinic Foundation

Donna M. Meyer, Ph.D.,
Immediate Past President
CHRISTUS Health

John V. Bergen, Ph.D.,
Executive Vice President

BOARD OF DIRECTORS

Susan Blonshine, RRT, RPFT, FAARC
TechEd

Kurt H. Davis, FCSMLS, CAE
Canadian Society for Medical Laboratory Science

Mary Lou Gantzer, Ph.D.
Dade Behring Inc.

Lillian J. Gill, M.S.
FDA Center for Devices and Radiological Health

Carolyn D. Jones, J.D., M.P.H.
AdvaMed

J. Stephen Kroger, M.D., MACP
COLA

Willie E. May, Ph.D.
National Institute of Standards and Technology

Gary L. Myers, Ph.D.
Centers for Disease Control and Prevention

Klaus E. Stinshoff, Dr.rer.nat.
Digene (Switzerland) Sàrl

James A. Thomas
ASTM International

Kiyoaki Watanabe, M.D.
Keio University School of Medicine

Judith A. Yost, M.A., M.T.(ASCP)
Centers for Medicare & Medicaid Services

Explore the Latest Offerings From CLSI!

As we continue to set the global standard for quality in laboratory testing, we are adding products and programs to bring even more value to our members and customers.



By becoming a CLSI member, your laboratory will join 1,600+ other influential organizations all working together to further CLSI's efforts to improve health care outcomes. You can play an active role in raising global laboratory testing standards—in your laboratory, and around the world.

Find out which membership option is best for you at www.clsi.org/membership.



Find what your laboratory needs to succeed! CLSI U provides convenient, cost-effective continuing education and training resources to help you advance your professional development. We have a variety of easy-to-use, online educational resources that make eLearning stress-free and convenient for you and your staff.

See our current educational offerings at www.clsi.org/education.



When laboratory testing quality is critical, standards are needed and there is no time to waste. eCLIPSE™ Ultimate Access, our cloud-based online portal of the complete library of CLSI standards, makes it easy to quickly find the CLSI resources you need.

Learn more and purchase eCLIPSE at cls.org/eCLIPSE.

For more information, visit www.clsi.org today.



950 West Valley Road, Suite 2500, Wayne, PA 19087 USA

P: +1.610.688.0100 **Toll Free (US):** 877.447.1888 F: +1.610.688.0700

E: customerservice@clsi.org www.clsi.org

ISBN 1-56238-550-X